
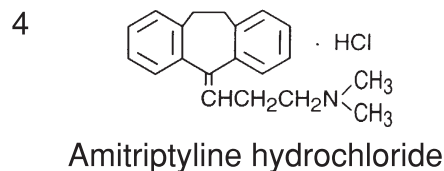
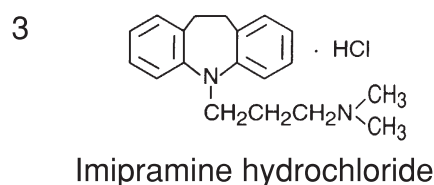
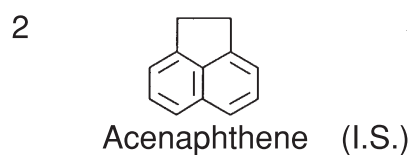
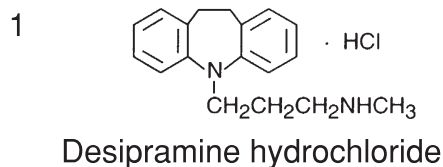


三環系抗うつ薬
Tricyclic antidepressants

A961015G

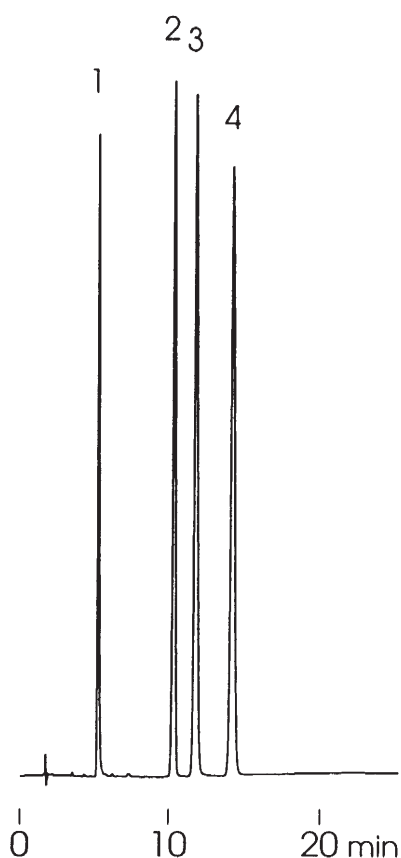

Column	: YMC-Pack Pro C18 (5 μm, 12nm) 150 × 4.6 mm I.D.
Eluent	: 20mM KH ₂ PO ₄ -K ₂ HPO ₄ (pH 6.9)/methanol (25/75)
Flow rate	: 1.0 mL/min
Temperature	: 37 °C
Detection	: UV at 254 nm, 0.16 AUFS
Injection	: 5 μL (0.1~0.5 mg/mL)