
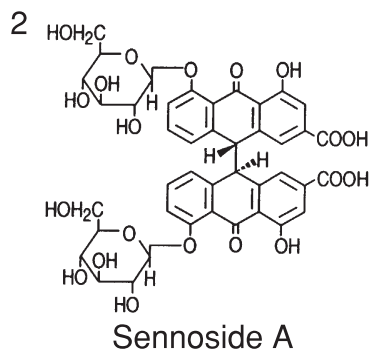
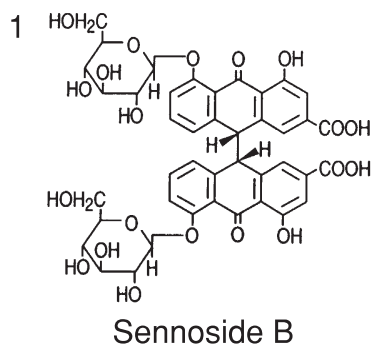


センノシド  
Sennosides

A961126B

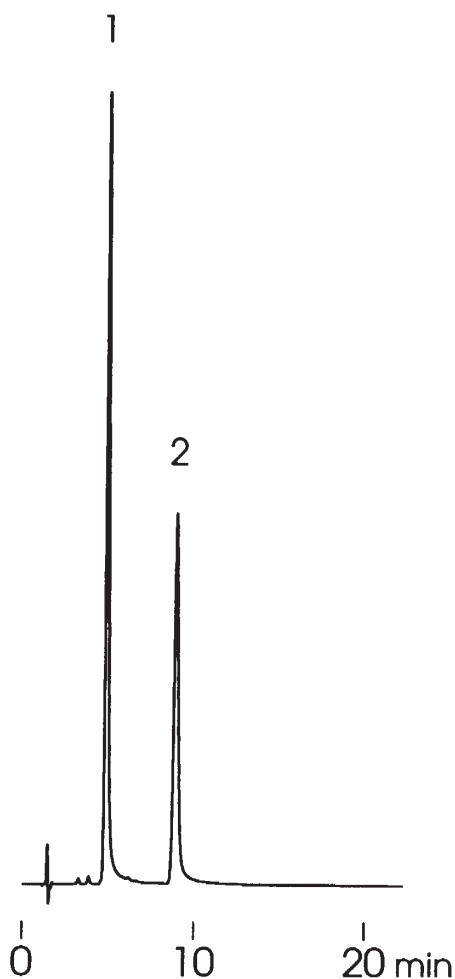

|             |   |
|-------------|---|
| Column      | : J'sphere ODS-H80 (4 $\mu$ m, 8nm)<br>150 $\times$ 4.6 mm I.D. |
| Eluent      | : methanol/water/acetic acid (30/65/5)                          |
| Flow rate   | : 1.0 mL/min  |
| Temperature | : 37 $^{\circ}$ C   |
| Detection   | : UV at 340 nm, 0.16 AUFS                                       |
| Injection   | : 10 $\mu$ L (0.1 mg/mL)  |