



ビタミンK
Vitamin K


A970728B

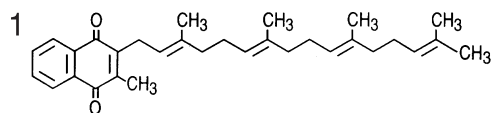

Menaquinone (Vitamin K2)

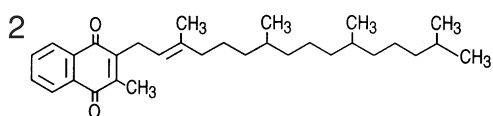

Phylloquinone (Vitamin K1)

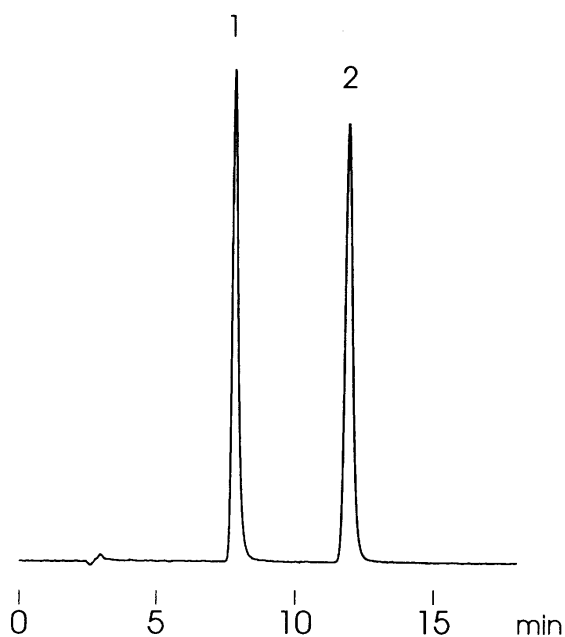

Column : YMC-Pack Pro C18 (5 μ m, 12nm)
150 \times 4.6 mm I.D.
Eluent : methanol/ethanol (5/1)
Flow rate : 0.8 mL/min
Temperature : 40 $^{\circ}$ C
Detection : FLS at Ex 320 nm, Em 430 nm
reduction column (platinum black, 10 x 4.0 mmI.D.)
Injection : 20 μ L (100 ng/mL)