
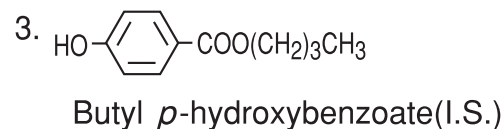
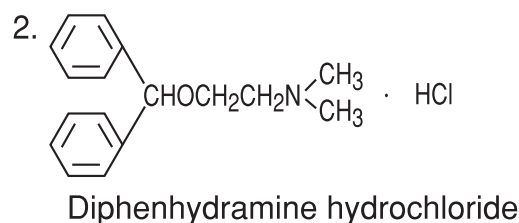
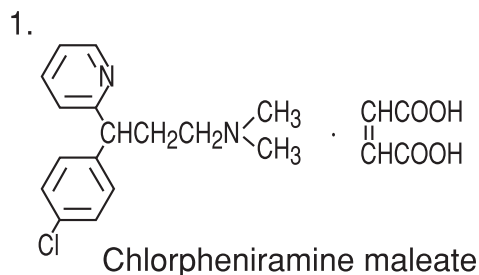


抗ヒスタミン薬  
Histamine antagonists

A980723A

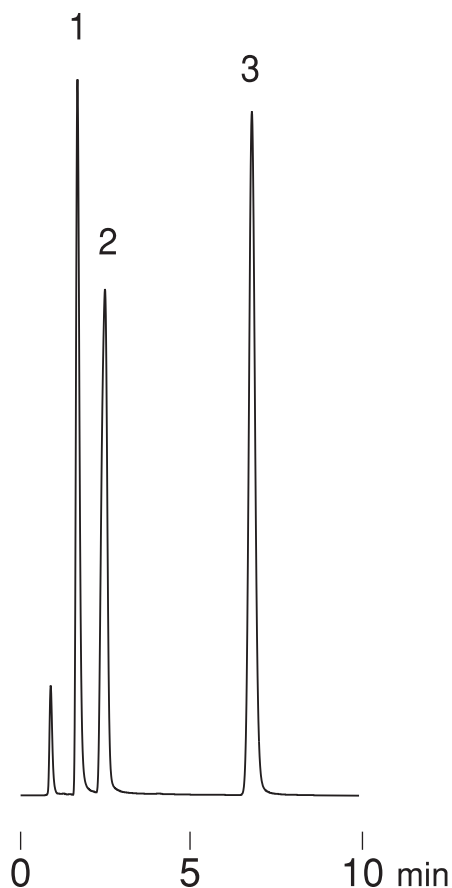

Column	: YMC-Pack Pro C18 (3 $\mu$ m, 12nm) 75 $\times$ 2.0 mm I.D.
Eluent	: acetonitrile / 50mM CH <sub>3</sub> COONH <sub>4</sub> (40/60)
Flow rate	: 0.2 mL/min
Temperature	: 37°C
Detection	: UV at 260 nm, 0.16 AUFS
Injection	: 2 $\mu$ L (0.05~2.0 mg/mL)