
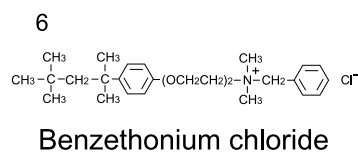
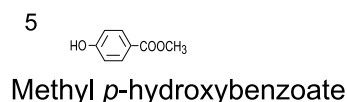
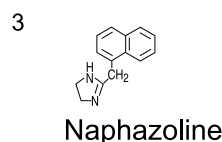
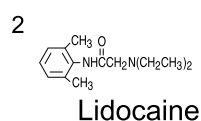



市販点鼻薬  
A commercial nasal spray

F040119B

A) Standard

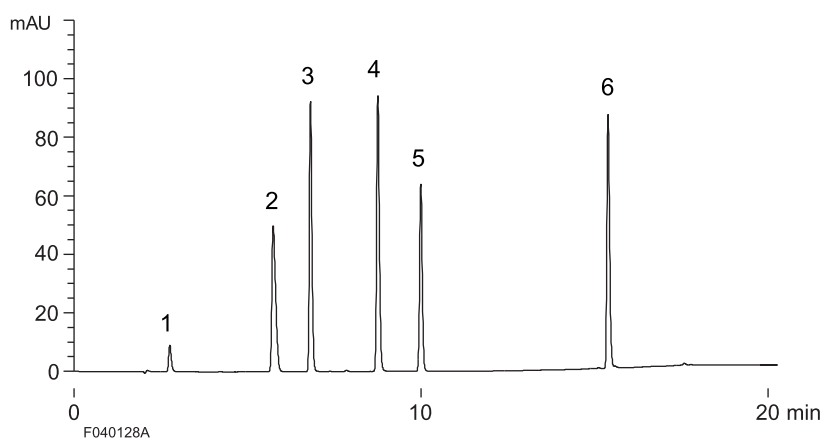

B) Nasal spray

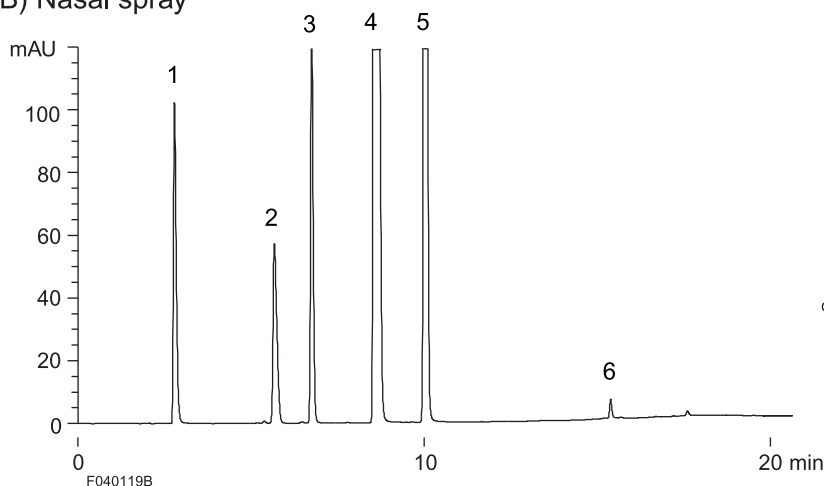

Column : YMC-Pack Pro C18 (5 μm, 12nm)  
150 × 4.6 mm I.D.

Eluent : A) 20mM KH<sub>2</sub>PO<sub>4</sub> - H<sub>3</sub>PO<sub>4</sub> (pH 2.5)  
B) methanol  
20-90%B(0-15min), 90%B(15-20min)

Flow rate : 1.0 mL/min

Temperature : 37 °C

Detection : UV at 260 nm

Injection : A) 5 μL (0.1~0.7 mg/mL)  
B) 3 μL (2 times dilution of a commercially available nasal spray with water)