


タンパク分子量マーカー  
Proteins for molecular weight marker

G911214C

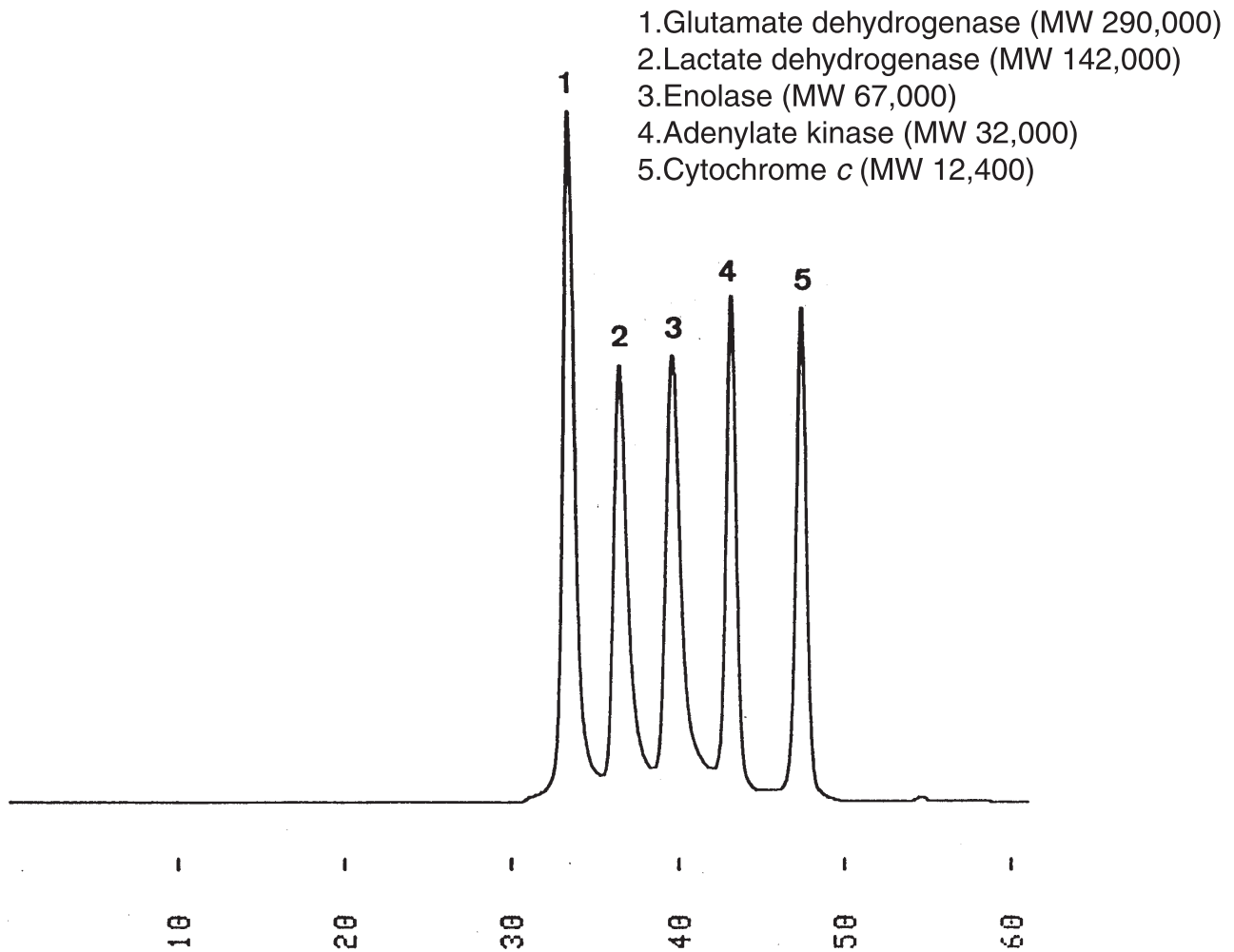

Column : YMC-Pack Diol-120 + Diol-300  
500 × 8.0 mm I.D. × 2

Eluent : 0.1M KH<sub>2</sub>PO<sub>4</sub>-K<sub>2</sub>HPO<sub>4</sub> (pH 7.0) containing 0.2M NaCl

Flow rate : 0.7 mL/min

Temperature : ambient (26 °C)

Detection : UV at 280 nm, 0.08 AUFS

Injection : 15 μL (100 μL/1 vial)

Sample : MW-Marker(HPLC)  
manufactured by ORIENTAL YEAST CO.,LTD.