



抗ヒスタミン薬  
Histamine antagonists

J010619A

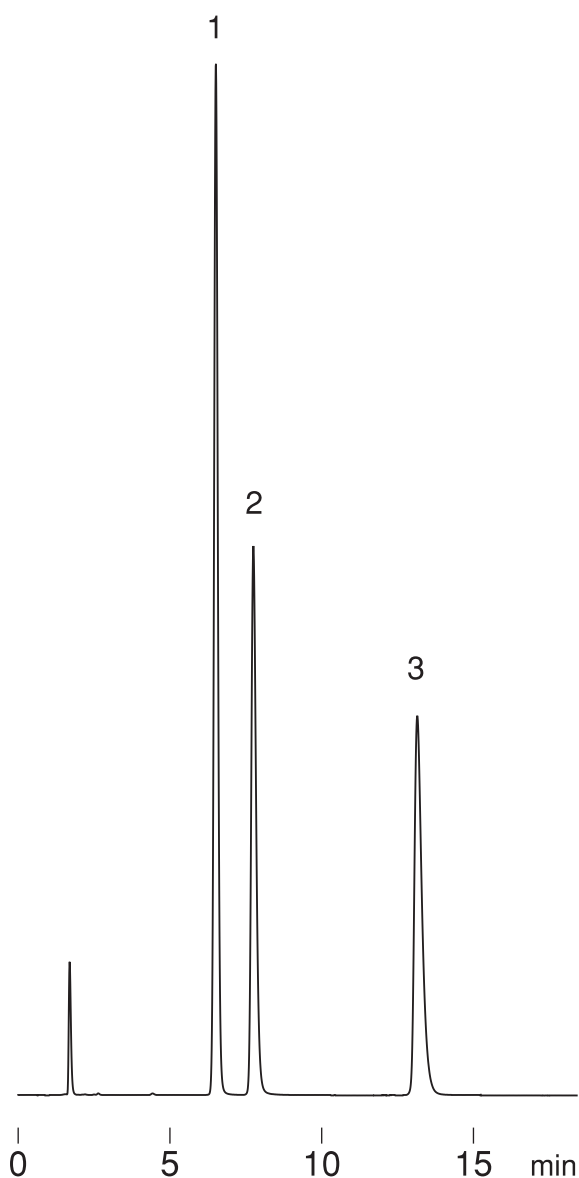


1.

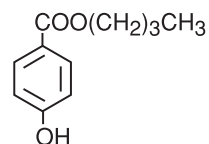


Butyl *p*-hydroxybenzoate (I.S.)

2.

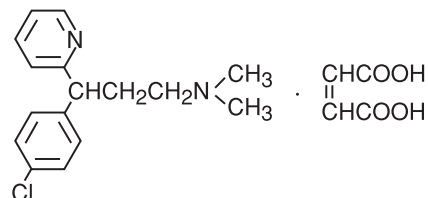

Chlorpheniramine maleate

3.

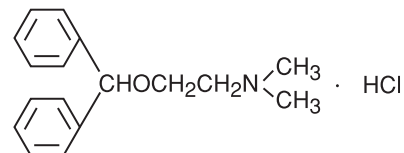

Diphenhydramine hydrochloride

Column	: YMC-Pack Pro C18 (5 μm, 12nm) 150 × 4.6 mm I.D.
Eluent	: 20mM KH <sub>2</sub> PO <sub>4</sub> -K <sub>2</sub> HPO <sub>4</sub> (pH 6.9) / methanol (35/65)
Flow rate	: 1.0 mL/min
Temperature	: 37 °C
Detection	: UV at 260 nm, 0.13 AUFS
Injection	: 5 μL (0.05~2.0 mg/mL)