


タンパク質
Proteins

N071016E

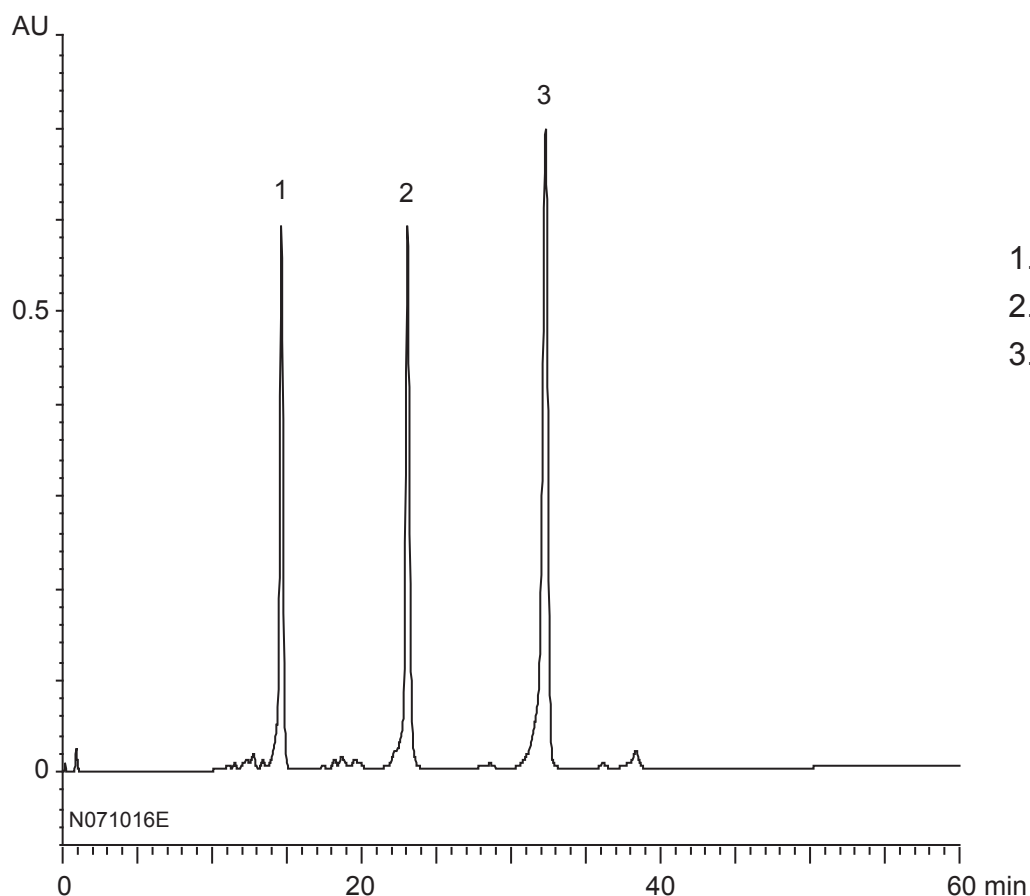

1. Ribonuclease A
2. Cytochrome c
3. Lysozyme

Column : BioPro IEX SP (5 μ m, 100 nm)
50 X 4.6 mm I.D.

Eluent : A) 20 mM KH_2PO_4 - K_2HPO_4 (pH 6.8)
B) 20 mM KH_2PO_4 - K_2HPO_4 (pH 6.8) containing 0.5 M NaCl
0-100%B (0-60 min)

Flow rate : 0.5 mL/min

Temperature : 25°C

Detection : UV at 220 nm

Injection : 20 μ L (0.5 mg/mL)