
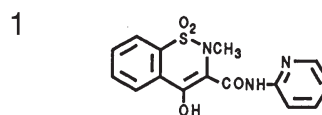



酸性非ステロイド性抗炎症薬

Acidic non-steroidal anti-inflammatory drugs


O950124A

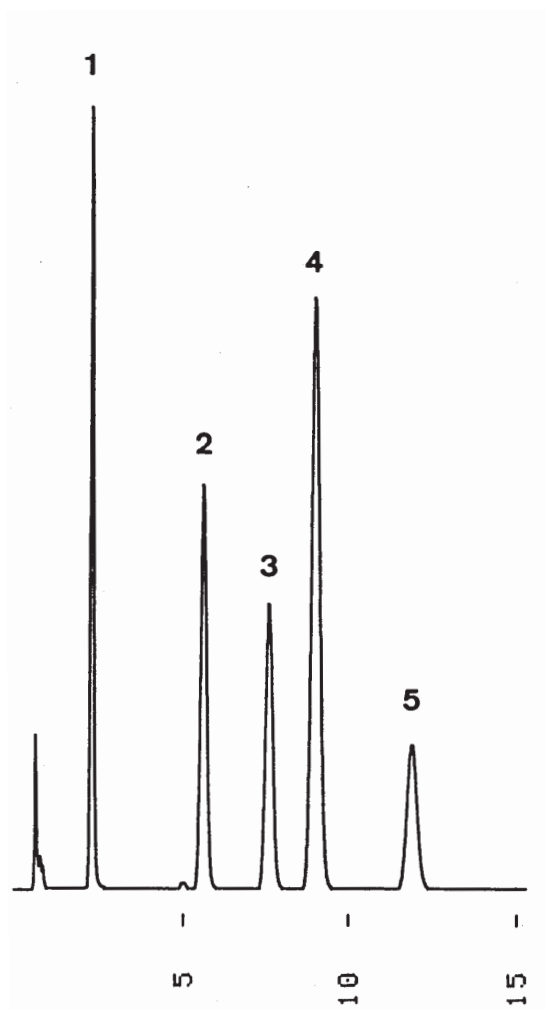


Piroxicam

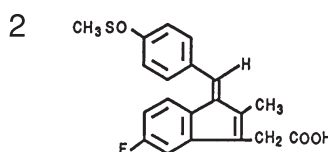


Sulindac

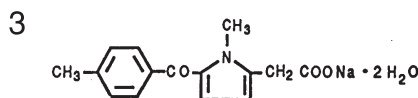

Tolmetin sodium dihydrate

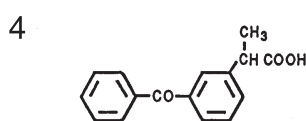

Ketoprofen

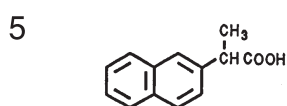

Naproxen

Column : J'sphere ODS-H80 (4 μ m, 8nm)
75 \times 4.6 mm I.D.
Eluent : 2-propanol/water/acetic acid (25/75/1)
Flow rate : 1.0 mL/min
Temperature : 37 $^{\circ}$ C
Detection : UV at 260 nm, 0.04 AUFS
Injection : 5 μ L (20 μ g/mL)