
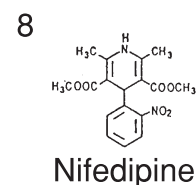
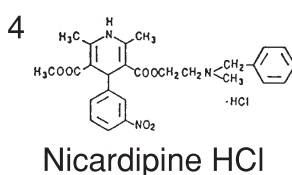
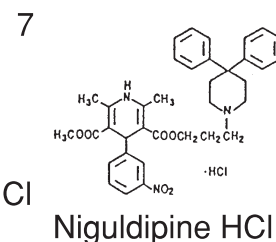
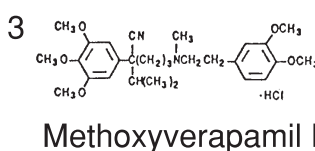
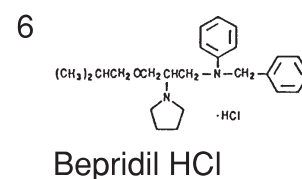
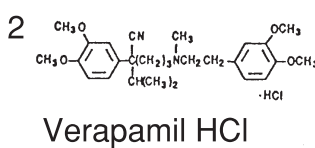
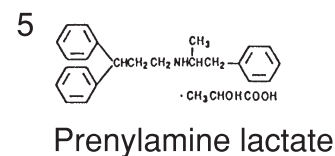
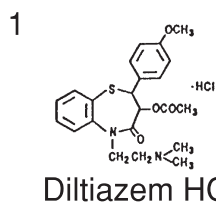


カルシウム拮抗薬  
Calcium antagonists

S930702A

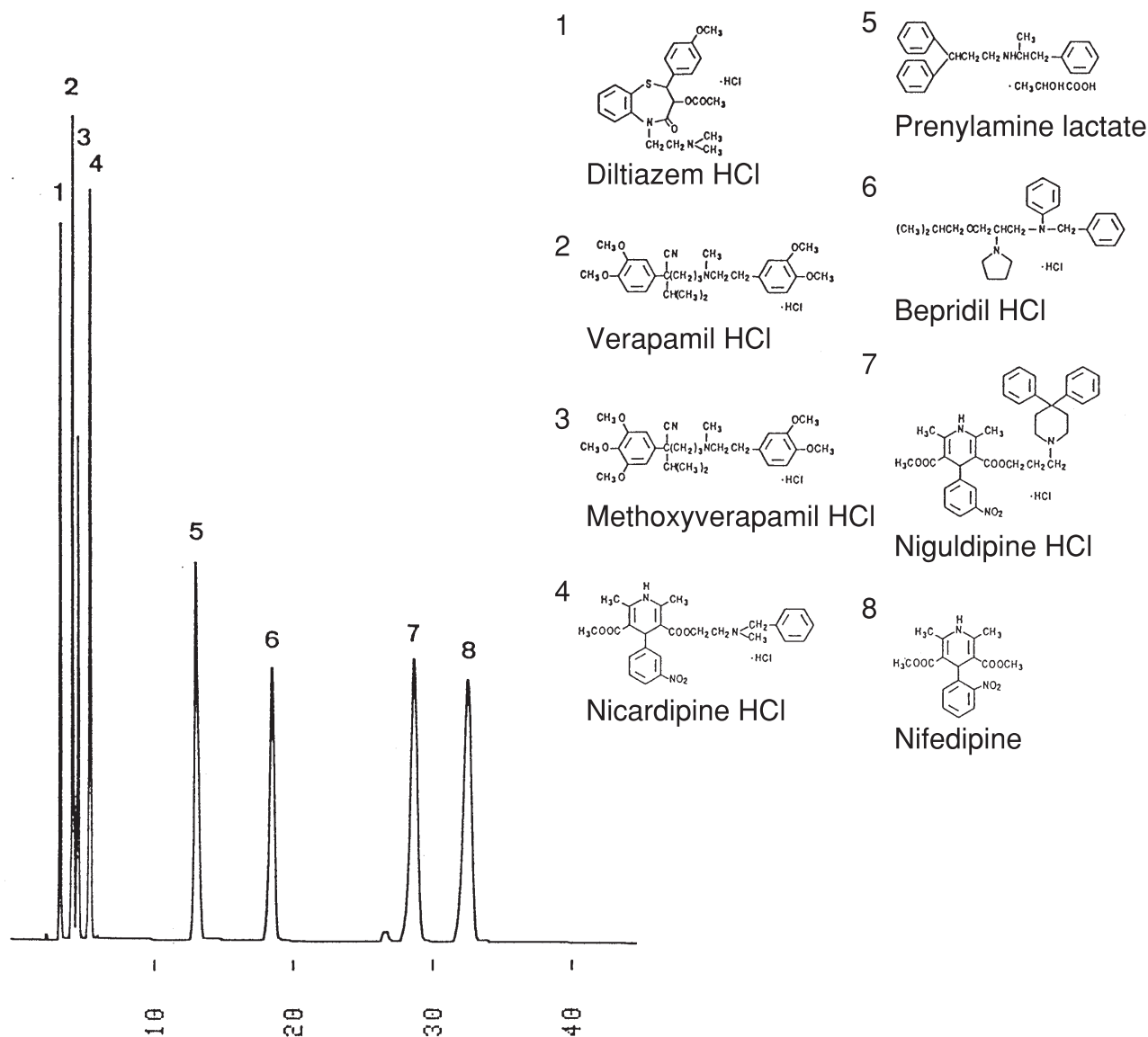

Column : J'sphere ODS-L80 (4  $\mu$  m, 8nm)  
150  $\times$  4.6 mm I.D.

Eluent : 20mM KH<sub>2</sub>PO<sub>4</sub> - H<sub>3</sub>PO<sub>4</sub> (pH 3.0)/acetonitrile/THF/TBA·Br  
(70/20/10/0.35, v/v/v/w)

Flow rate : 0.8 mL/min

Temperature : 30  $^{\circ}$ C

Detection : UV at 270 nm, 0.064 AUFS

Injection : 8  $\mu$ L (0.05~1.4 mg/mL)