



不飽和脂肪酸メチルエステル
Unsaturated fatty acid methyl esters


S931012C

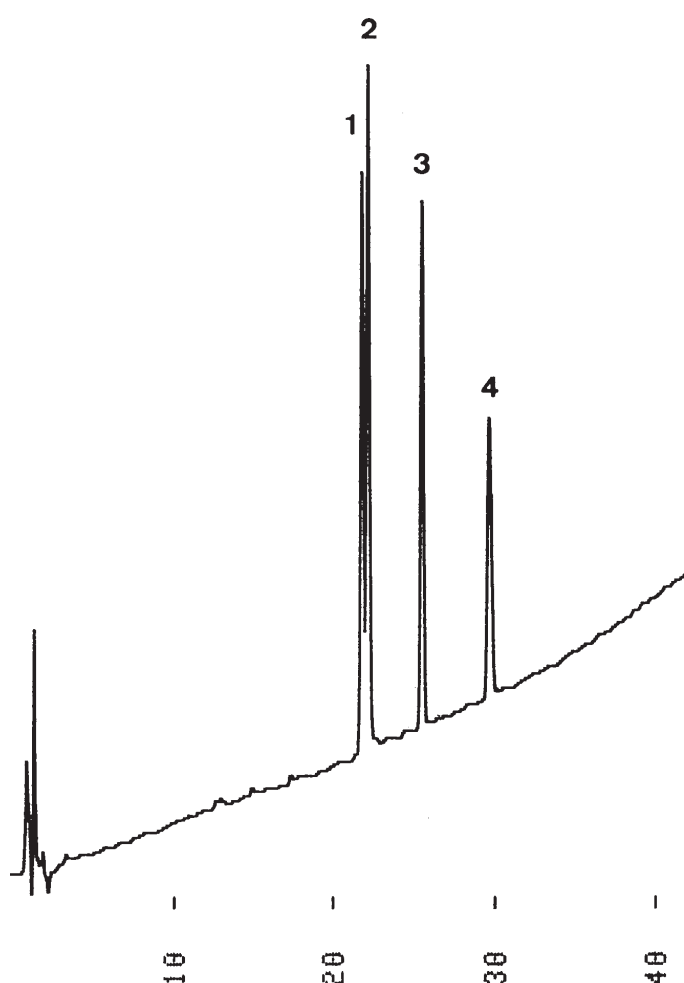

1. Methyl linolenate

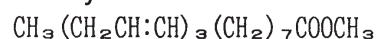

2. Methyl γ -linolenate

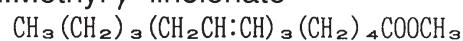

3. Methyl linoleate



4. Methyl oleate



Column : J'sphere ODS-H80 (4 μ m, 8nm)
150 \times 4.6 mm I.D.
Eluent : A)acetonitrile/THF/water (10/30/60)
B)acetonitrile/THF/water (65/30/5)
20-100%B (0-40min, linear)
Flow rate : 1.0 mL/min
Temperature : 37 $^{\circ}$ C
Detection : UV at 210 nm, 0.51 AUFS
Injection : 18 μ L (0.25~4 μ L/mL)