
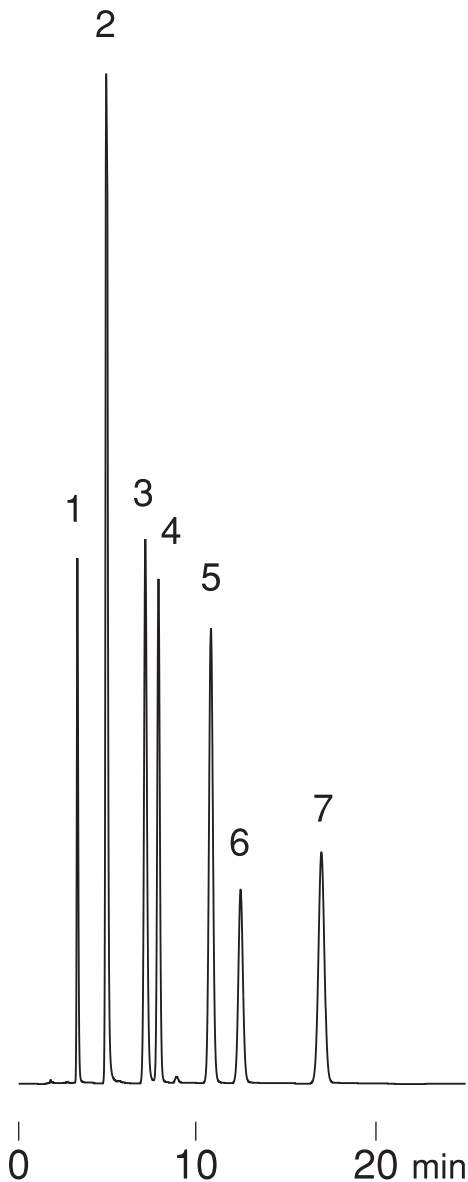
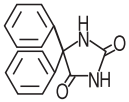
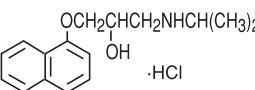
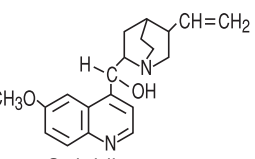
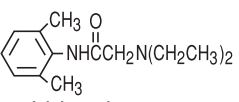
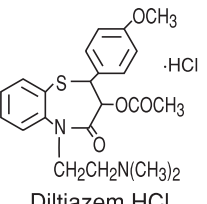
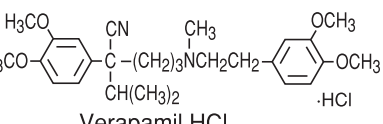


不整脈治療薬
Antiarrhythmic drugs

S990817C

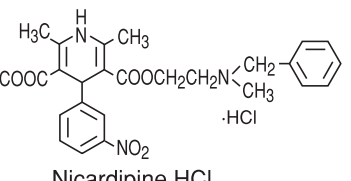

1. 
Phenytoin
2. 
Propranolol HCl
3. 
Quinidine
4. 
Lidocaine
5. 
Diltiazem HCl
6. 
Verapamil HCl
7. 
Nicardipine HCl

Column	: YMC-Pack Pro C4 (5 μ m, 12nm) 150 \times 4.6 mm I.D.
Eluent	: 20mM KH ₂ PO ₄ -K ₂ HPO ₄ (pH 6.9) / methanol (40/60)
Flow rate	: 1.0 mL/min
Temperature	: 37 $^{\circ}$ C
Detection	: UV at 220 nm, 0.26 AUFS
Injection	: 5 μ L (0.05~0.1 mg/mL)