
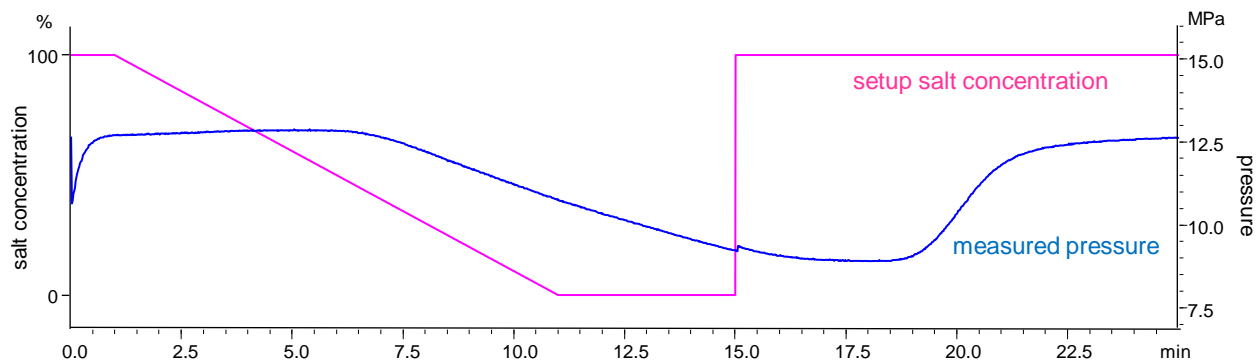


抗体医薬品
Antibody drugs

W191118F

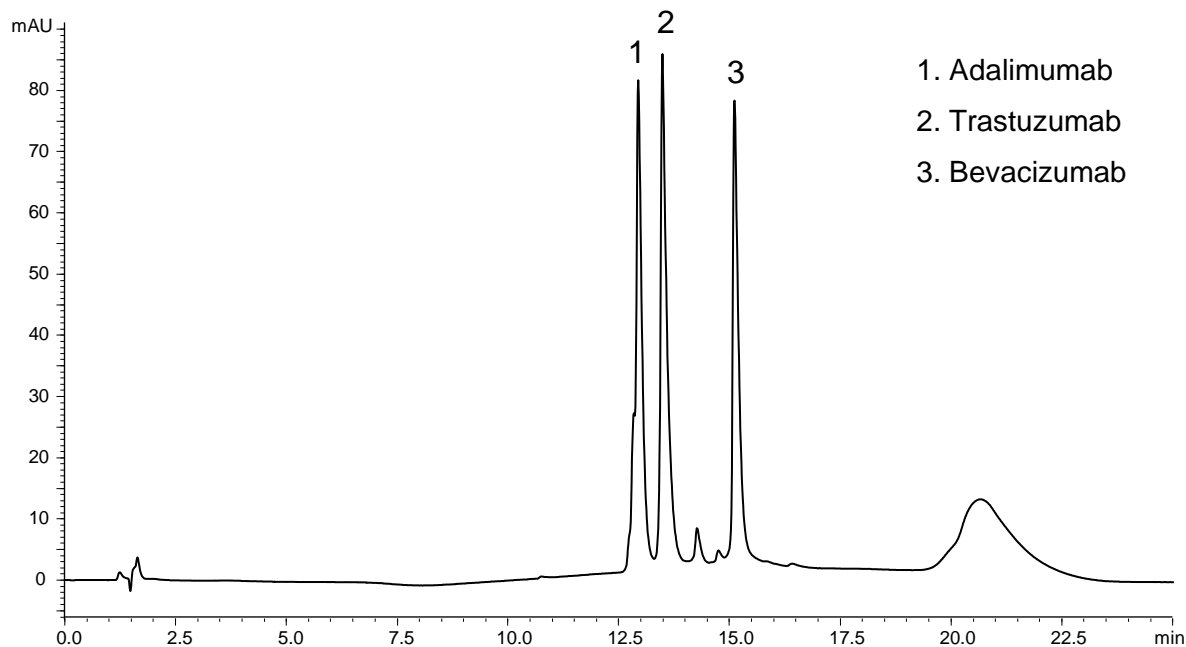

Column : BioPro HIC HT (2.3 μ m)

100 X 4.6 mm I.D.

Eluent : A) 100 mM NaH_2PO_4 - Na_2HPO_4 containing 2.0 M $(\text{NH}_4)_2\text{SO}_4$ (pH 7.0)

B) 100 mM NaH_2PO_4 - Na_2HPO_4 (pH 7.0)

0%B (0-1 min), 0-100%B (1-11 min), 100%B (11-15 min), 0%B (15-25 min)

Flow rate : 0.5 mL/min

Temperature : 25°C

Detection : UV at 280 nm

Injection : 15 μ L (0.5 mg/mL)