

Preparative HPLC systems

K-Prep FC

K-Prep FC is a preparative HPLC system that allows seamless scale-up from small scale production up to industrial scaling-up processes.



Features


- Smooth transition from laboratory scale to industrial scale
- Usable for production purpose
- All-in-one system equipped with HPLC pumps to load samples and fraction collectors
- Strong support for preparative purifications provided by fully-automatic operation by PC
- Compliance with IQ/OQ validation and CSV

Channel diagram

MARK	NAME
E	Inlet Port
F	Fraction Port
P	Pump
C	Column
3PBV	3way Pneumatic Ball Valve
4PBV	4way Pneumatic Ball Valve
4MBV	4way Manual Ball Valve
PE	Pressure Sensor
UV	UV Detector

MARK	NAME
TE	Temperature Sensor
PH	pH Sensor
EC	Conductivity Sensor
MIX	Mixer
RV	Relief Valve
BPV	Back Pressure Valve
PG	Pressure Gauge
DF	Disk Filter
WP	Plunger Washing Pump

* Gradient mode
 optional

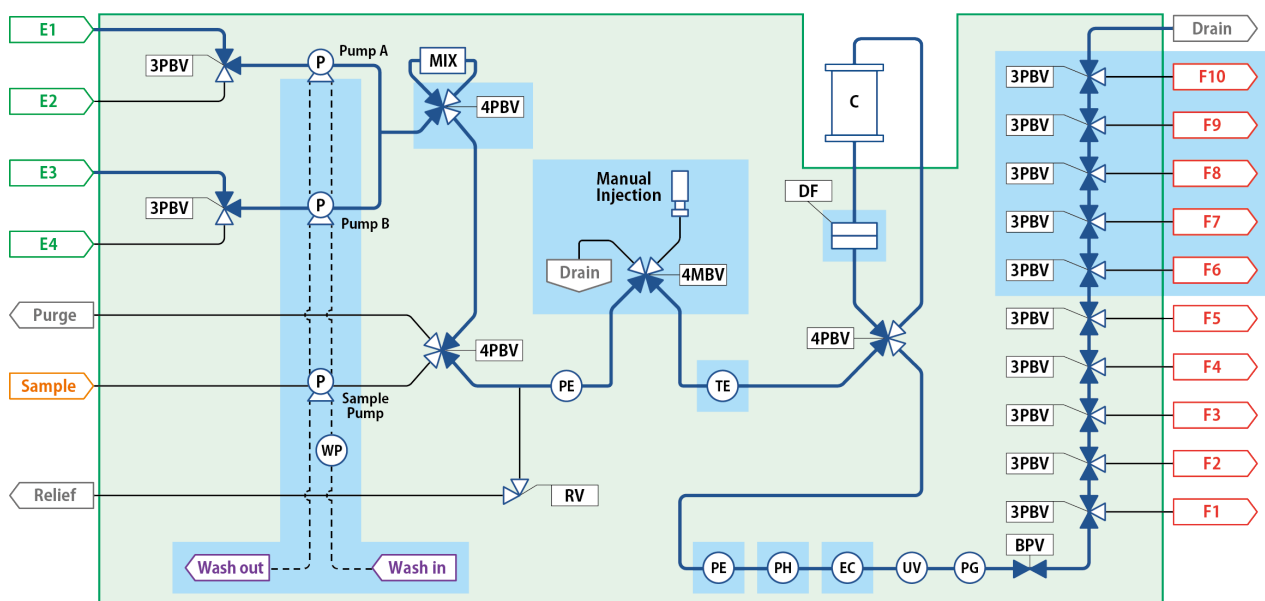

Specification

	FC750S	FC750G	FC1500S	FC1500G
Gradient capability	no	yes	no	yes
Flow rate range	1 - 750 mL/min		1 - 1500 mL/min	
Flow rate accuracy	±2% (75 mL/min and above)		±2% (150 mL/min and above)	
Maximum pressure	10 MPa			
Fraction port	1/4", SUS316, Swagelok fitting		1/4", SUS316, Swagelok fitting	
Inlet port	1/4", SUS316, Swagelok fitting		3/8", SUS316, Swagelok fitting	
Solvent channel	2ch	4ch (2 × 2ch)	2ch	4ch (2 × 2ch)
Fraction collector	5-channel switching valve method + drain			
Sample pump	1 - 300 mL/min			
Detector	UV (200 - 400 nm), 3 variable-wavelengths measurable			
Control and display	laptop PC (changeable FA PC) and programmable logic controller			
Utility	AC200V 50/60Hz			
	0.6 kW	1.1 kW	0.6 kW	1.1 kW
	0.5 MPa dry air			
Dimension	850 (W) × 970 (D) × 1300 (H) mm			

Before use (installation, operation, maintenance or check-up) of our products, instruction manuals should be carefully read and understood, and safety rules and precautions followed as outlined in the manuals.

Worldwide Availability

YMC America, Inc.
www.ymcamerica.com

YMC Europe GmbH
ymc.eu

YMC Switzerland LLC
www.ymc-schweiz.ch

YMC Shanghai Rep. Office
www.ymcchina.com

YMC India Pvt. Ltd.
www.ymcindia.com

YMC Korea Co., Ltd.
www.ymckorea.com

YMC Taiwan Co., Ltd.
www.ymctaiwan.com

YMC Singapore Tradelinks Pte. Ltd.
ymc.sg

YMC CO., LTD.

YMC Karasuma-Gojo Bldg., 284 Daigo-cho, Karasuma Nishiiru Gojo-dori, Shimogyo-ku, Kyoto, 600-8106, Japan
TEL: +81-75-342-4515 FAX: +81-75-342-4550
www.ymc.co.jp sales@ymc.co.jp