


メトプロロール


Metoprolol

H131008D

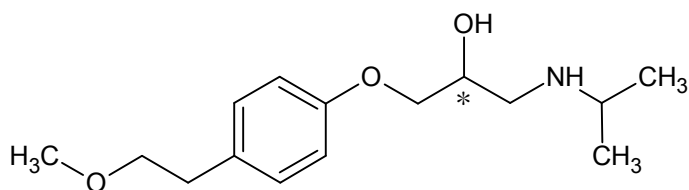

Metoprolol

|          |      |
|----------|------|
| $\alpha$ | 1.45 |
| $R_s$    | 5.81 |

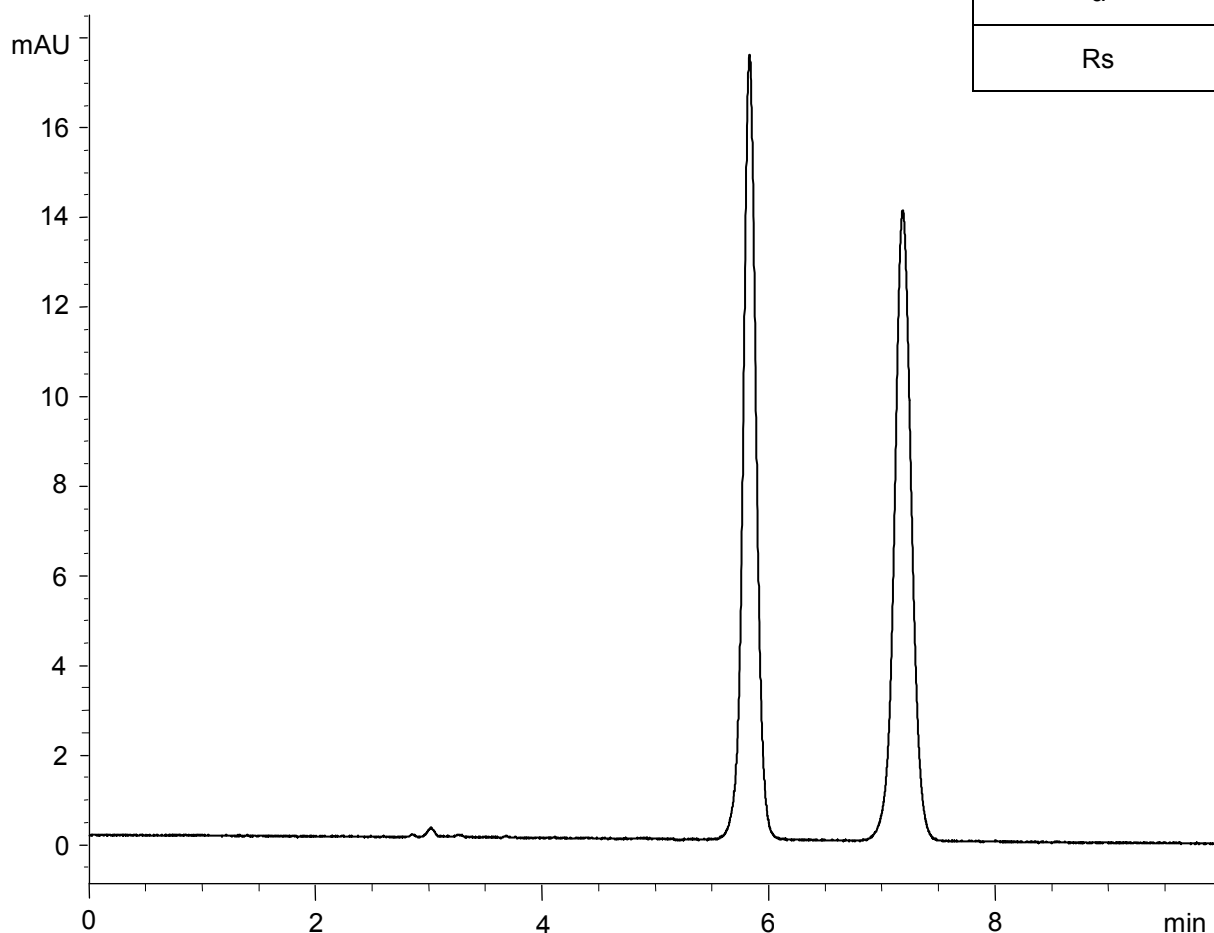


Column : YMC CHIRAL Amylose-C (5  $\mu$ m)  
250 X 4.6 mm I.D.  
Eluent : *n*-hexane/ethanol/ethanolamine (80/20/0.1)  
Flow rate : 1.0 mL/min  
Temperature : 25°C  
Detection : UV at 254 nm  
Injection : 5  $\mu$ L (1 mg/mL)


メトプロロール

Metoprolol

H131009B

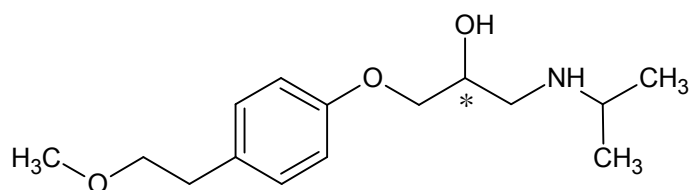

Metoprolol

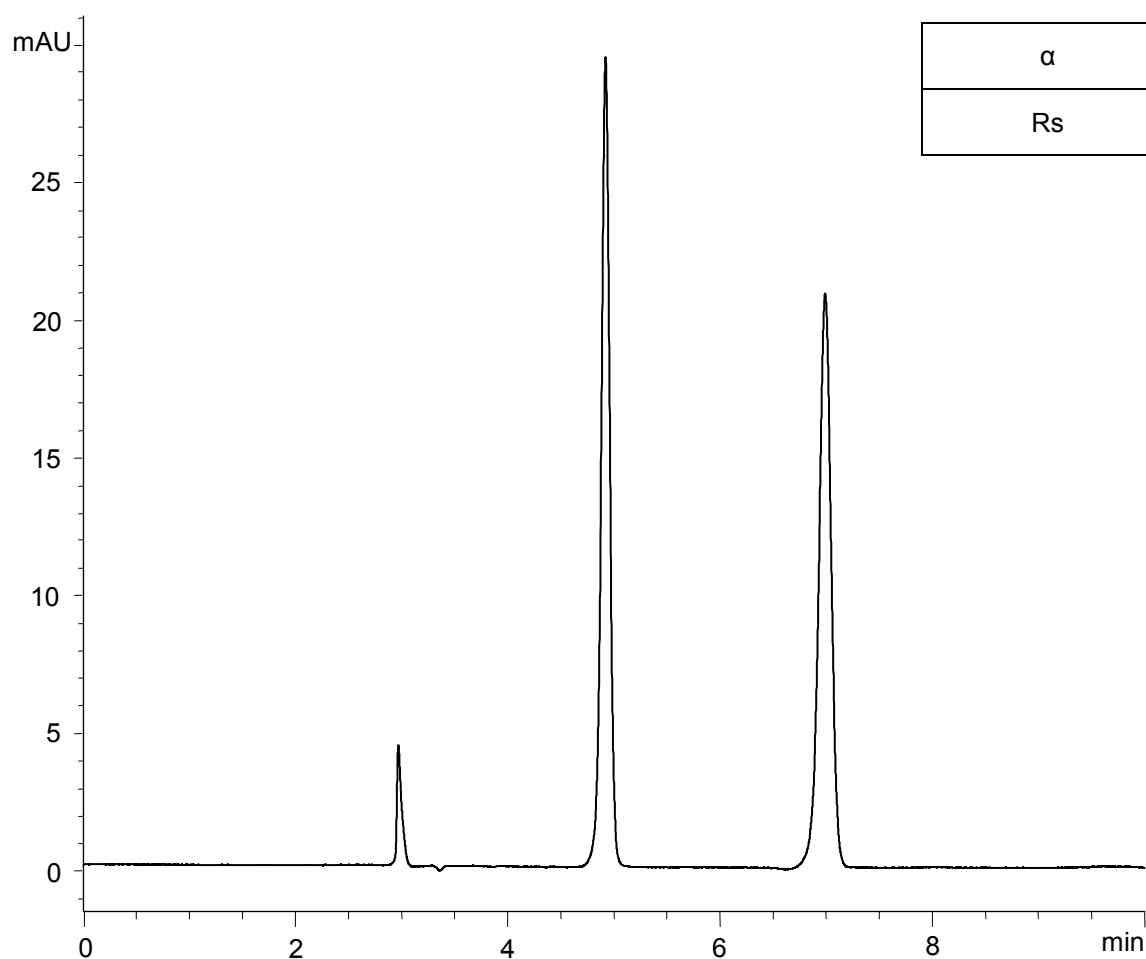

Column : YMC CHIRAL Cellulose-C (5  $\mu$ m)  
 250 X 4.6 mmI.D.  
 Eluent : *n*-hexane/ethanol/ethanolamine (90/10/0.1)  
 Flow rate : 1.0 mL/min  
 Temperature : 25°C  
 Detection : UV at 254 nm  
 Injection : 5  $\mu$ L (1 mg/mL)