


プラスミドpBR322とその制限酵素*Hae* III分解物

Plasmid pBR322 and pBR322 *Hae* III restriction fragments

P080617A

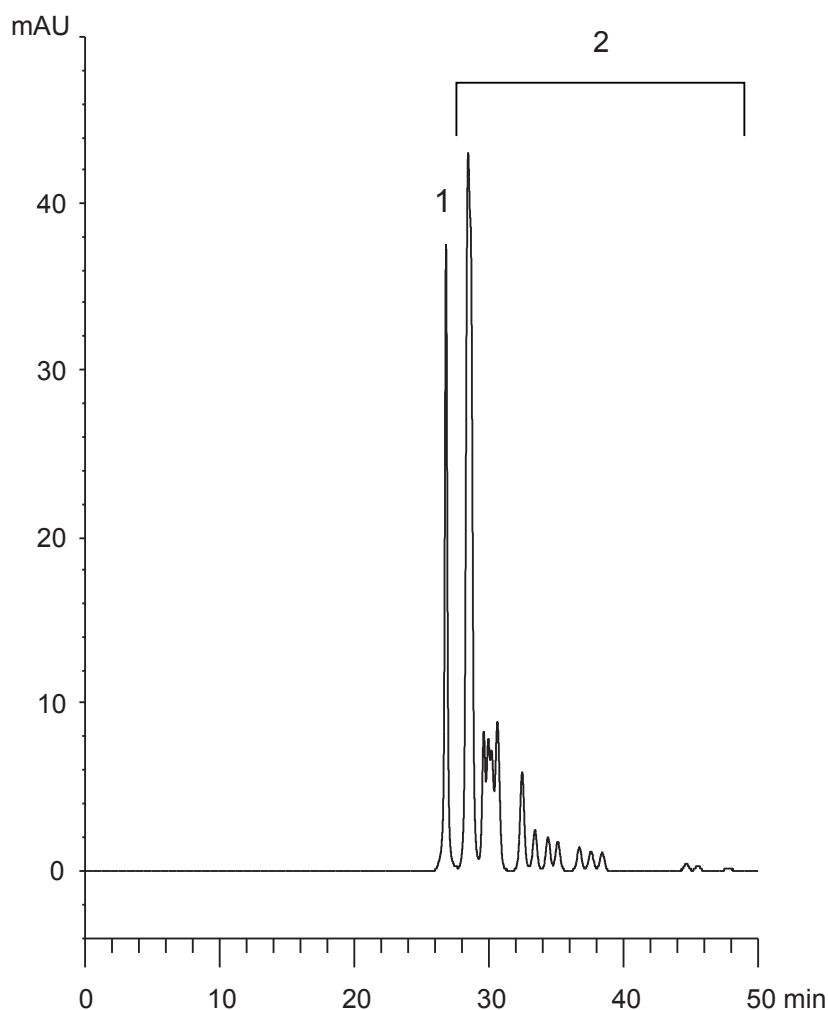

1. Plasmid pBR322 (4,361 bp)

2. Plasmid pBR322 *Hae* III digest (8-587 bp)

Column : YMC-Pack Diol-300 + Diol-200  
500 X 8.0 mmI.D. X 2

Eluent : 0.1 M KH<sub>2</sub>PO<sub>4</sub>-K<sub>2</sub>HPO<sub>4</sub> (pH 7.0) containing 0.2 M NaCl

Flow rate : 0.7 mL/min

Temperature : ambient (25°C)

Detection : UV at 260 nm

Injection : 10 μL

Sample : Plasmid pBR322 (0.03 mg/mL)  
Plasmid pBR322 *Hae* III digest (0.13 mg/mL)