


アミノ酸  
Amino acids

A000419B

1. L-Valine
2. L-Tyrosine
3. L-Phenylalanine
4. L-Tryptophan

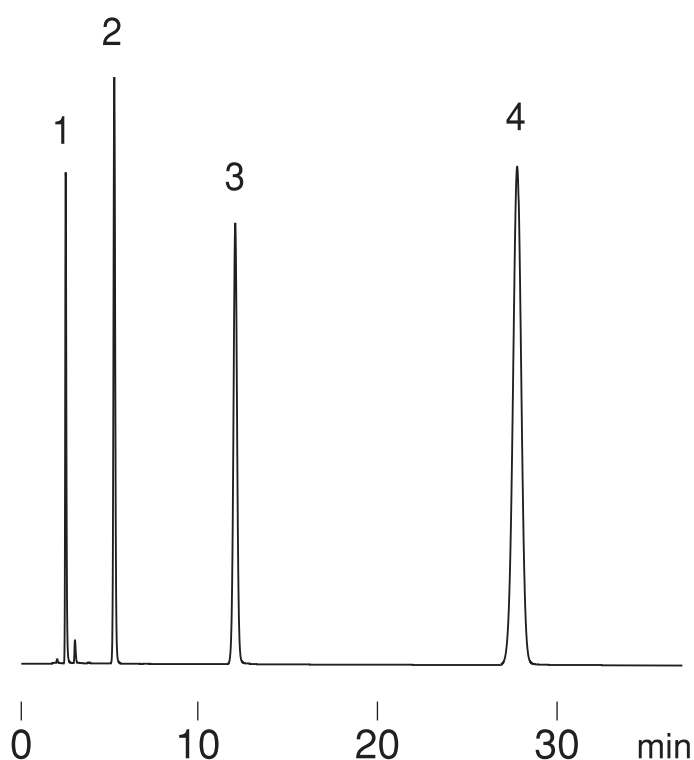

Column : Hydrosphere C18 (5  $\mu$ m, 12nm)  
150  $\times$  4.6 mm I.D.  
Eluent : 20mM KH<sub>2</sub>PO<sub>4</sub>  
Flow rate : 1.0 mL/min  
Temperature : 37  $^{\circ}$ C  
Detection : UV at 210 nm, 0.16 AUFS  
Injection : 5  $\mu$ L (0.05~1.6 mg/mL)