


有機酸
Organic acids

A001218A

1. Tartaric acid
2. L-Malic acid
3. Lactic acid
4. Acetic acid
5. Citric acid
6. Succinic acid

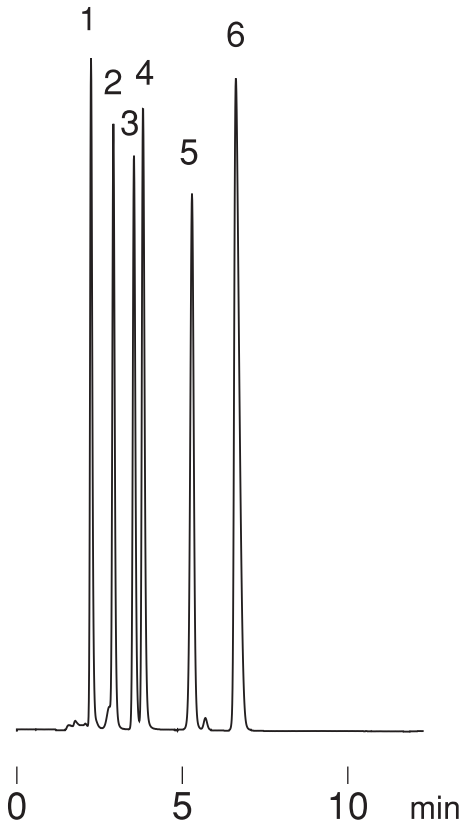

Column : Hydrosphere C18 (5 μ m, 12nm)
150 \times 4.6 mm I.D.
Eluent : 20mM NaH₂PO₄ - H₃PO₄ (pH 2.8)
Flow rate : 1.0 mL/min
Temperature : 30 $^{\circ}$ C
Detection : UV at 220 nm, 0.064 AUFS
Injection : 10 μ L (0.25~2.0 mg/mL)