



シタロプラム臭化水素酸塩  
Citalopram hydrobromide

AB201015A

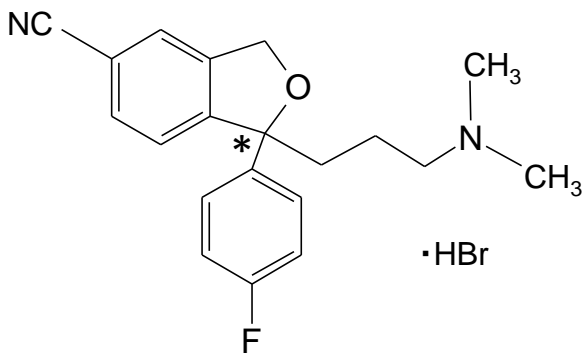

|          |     |
|----------|-----|
| $\alpha$ | 1.4 |
| $R_s$    | 6.0 |

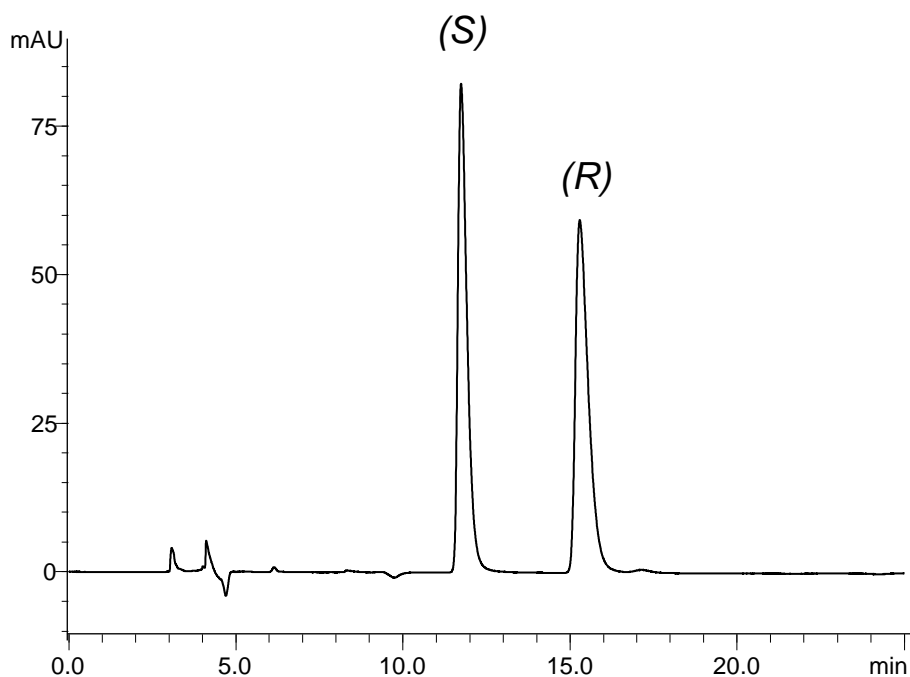

Column : CHIRAL ART Cellulose-SZ (5  $\mu$ m)  
250 X 4.6 mm I.D.  
Eluent : *n*-hexane/2-propanol/diethylamine (90/10/0.1)  
Flow rate : 1.0 mL/min  
Temperature : 25°C  
Detection : UV at 270 nm  
Injection : 7  $\mu$ L (1.0 mg/mL)