


トリペプチド  
Tripeptides

AB211108A

1. Val-Pro-Pro  
2. Ile-Pro-Pro

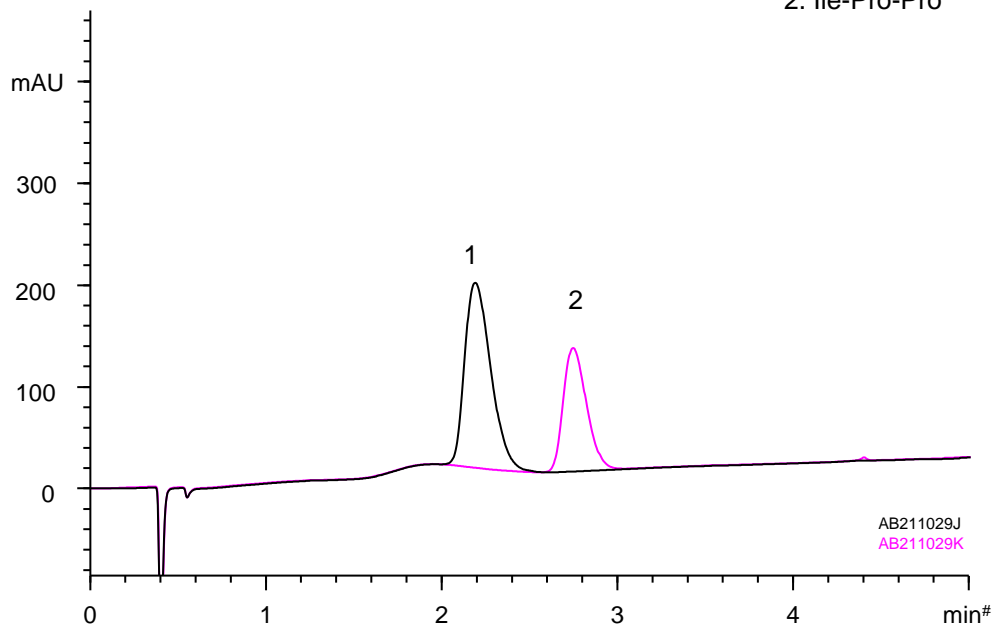

AB211029J  
AB211029K

Column	: YMC-Triart C18 (1.9 $\mu$ m, 12 nm) 50 X 2.1 mmI.D.
Eluent	: A) water/formic acid (100/0.1) B) acetonitrile/formic acid (100/0.08) 0-20%B (0-5 min)
Flow rate	: 0.4 mL/min
Temperature	: 60°C
Detection	: UV at 210 nm
Injection	: 2 $\mu$ L (0.1 mg/mL)