


配位性化合物
Coordination compounds


AS231120A

Accura Triart C18 (bioinert hardware)

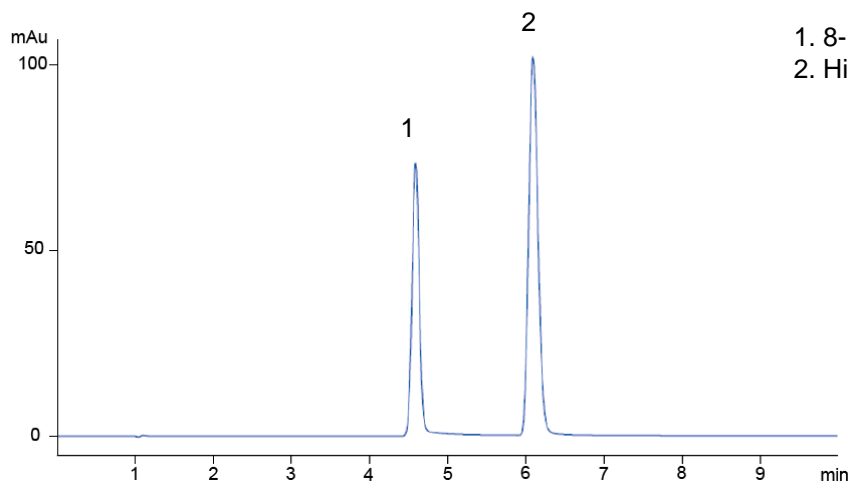

1. 8-Hydroxyquinoline (8-Quinololinol)
2. Hinokitiol

YMC-Triart C18 (stainless-steel hardware)

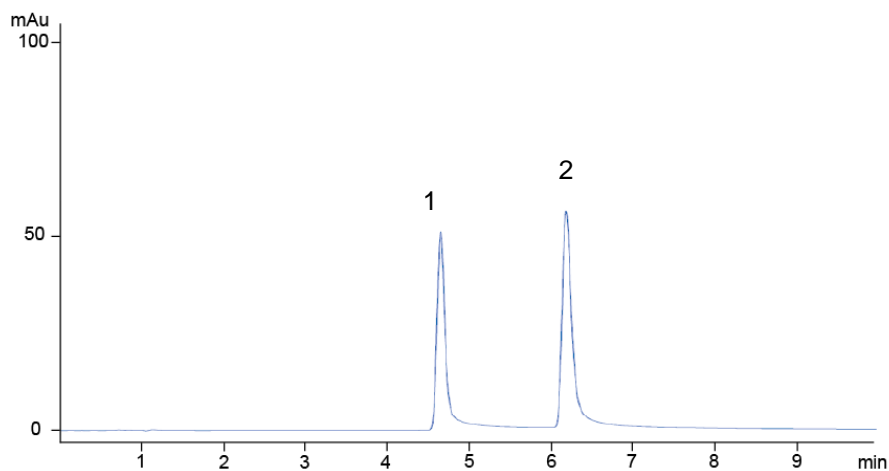

Courtesy of YMC America, Inc.

Column	: 3 μ m, 12 nm 100 X 4.6 mm I.D.
Eluent	: acetonitrile/water (40/60)
Flow rate	: 1.0 mL/min
Temperature	: 30°C
Detection	: UV at 310 nm
Injection	: 10 μ L (0.1 mg/mL)