


イマザリル  
Imazalil

B931007A

1.Imazalil

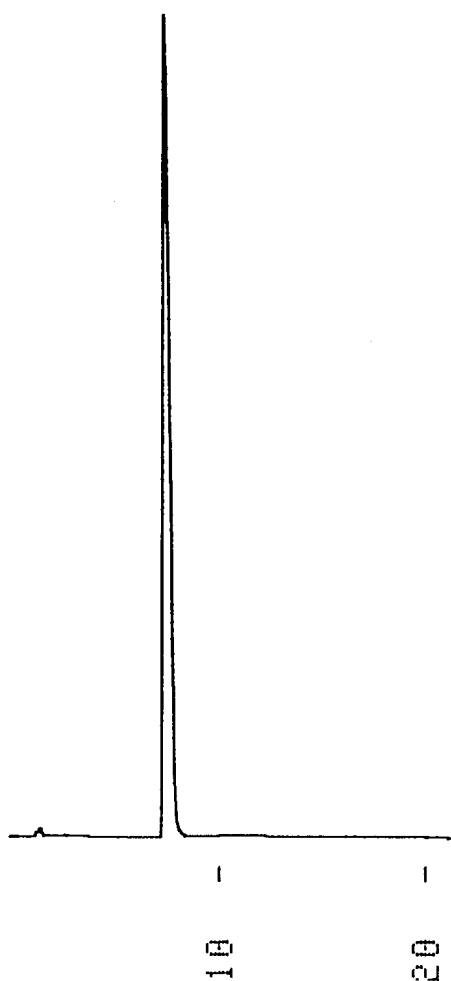

Column : J'sphere ODS-M80 (4  $\mu$  m, 8nm)  
150  $\times$  4.6 mm I.D.  
Eluent : acetonitrile/20mM NH<sub>4</sub>H<sub>2</sub>PO<sub>4</sub> (50/50)  
Flow rate : 1.0 mL/min  
Temperature : 37  $^{\circ}$ C  
Detection : UV at 230 nm, 0.04 AUFS  
Injection : 10  $\mu$ L (0.05 mg/mL)