
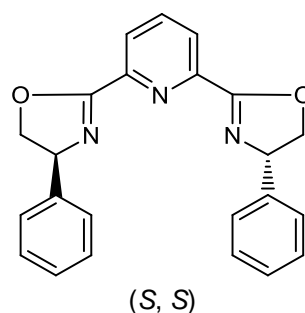
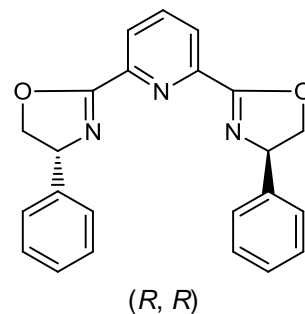


2,6-ビス(4-フェニル-2-オキサゾリン-2-イル)ピリジン
2,6-Bis(4-phenyl-2-oxazolin-2-yl)pyridine

E170801B

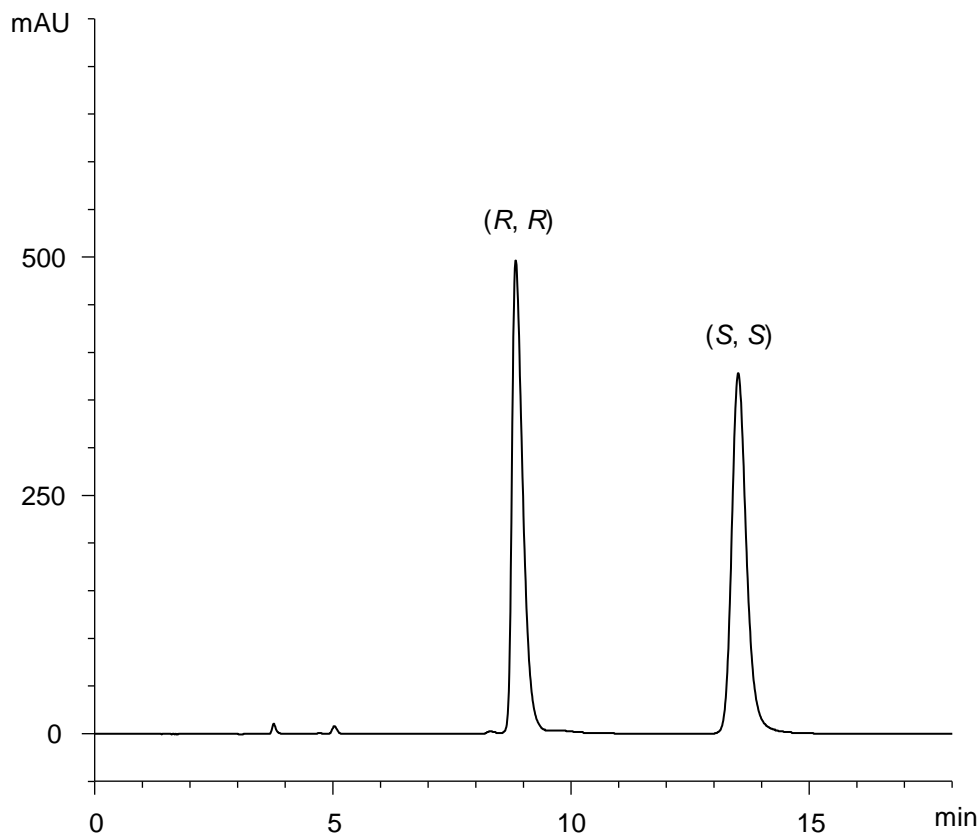

α	1.7
R_s	9.5

Column : CHIRAL ART Cellulose-C (5 μ m)
250 X 4.6 mm I.D.
Eluent : *n*-hexane/ethanol (40/60)
Flow rate : 1.0 mL/min
Temperature : 25°C
Detection : UV at 285 nm
Injection : 10 μ L (0.5 mg/mL)
Sample : Supplied by Tokyo Chemical Industry Co., Ltd.