
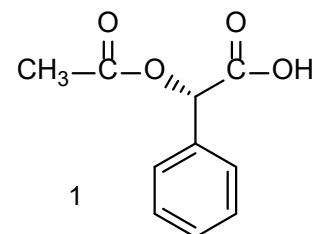



O-アセチルマンデル酸

O-Acetylmandelic Acid

E170904B

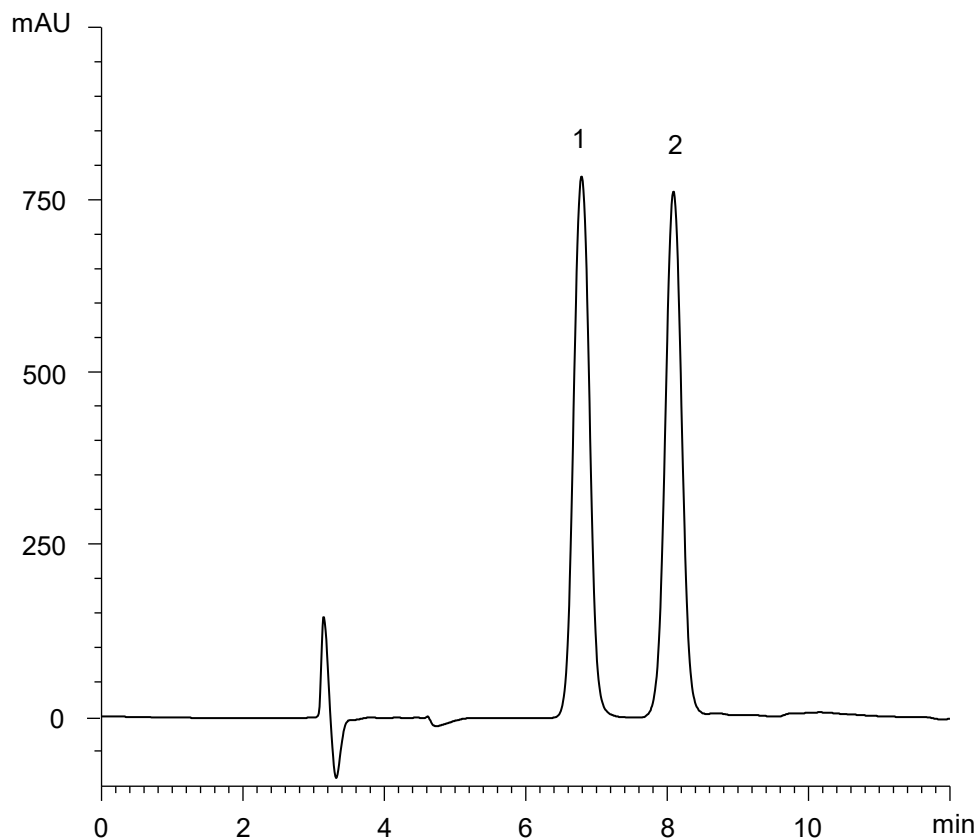

(+)-O-Acetyl-L-mandelic Acid

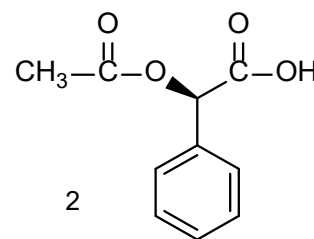

(-)-O-Acetyl-D-mandelic Acid

$\alpha$	1.3
$R_s$	3.0

Column : CHIRALART Cellulose-SC (5  $\mu$ m)  
250 X 4.6 mm I.D.

Eluent : *n*-hexane/2-propanol/trifluoroacetic acid (90/10/0.1)

Flow rate : 1.0 mL/min

Temperature : 25°C

Detection : UV at 210 nm

Injection : 10  $\mu$ L (0.5 mg/mL)

Sample : Supplied by Tokyo Chemical Industry Co., Ltd.