


アンジオテンシン
Angiotensins

G910520A

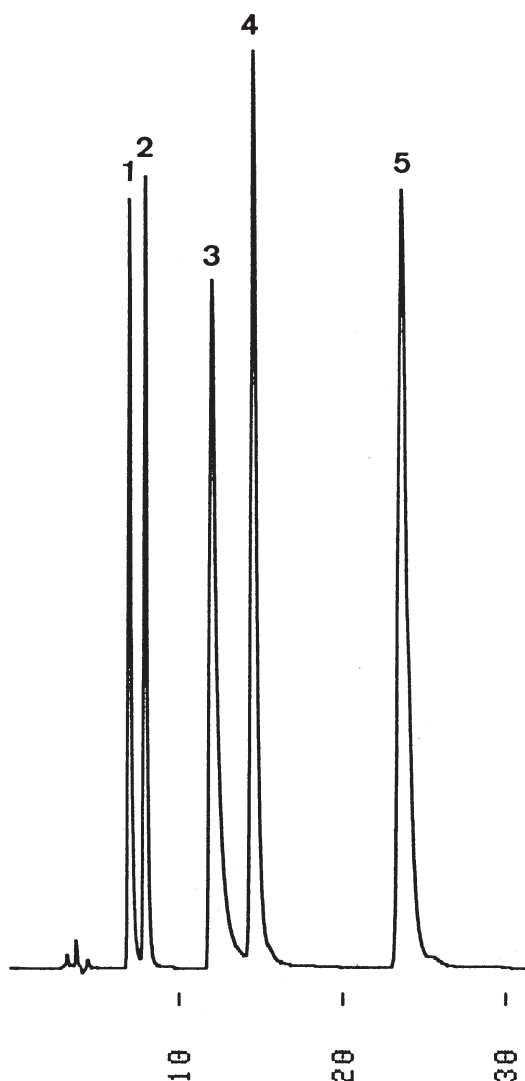

1. Angiotensin III
2. Angiotensin II
3. [Asn1, Val5] Angiotensin I
4. [Val5] Angiotensin I
5. Angiotensin I

Column : YMC-Pack ODS-AM (5 μ m, 12nm)
150 \times 4.6 mm I.D.
Eluent : acetonitrile/20mM KH₂PO₄ (20/80)
Flow rate : 0.5 mL/min
Temperature : 30 $^{\circ}$ C
Detection : UV at 220 nm, 0.08 AUFS
Injection : 8 μ L (0.05~0.15 mg/mL)