


カゼイン  
Casein

G920123B

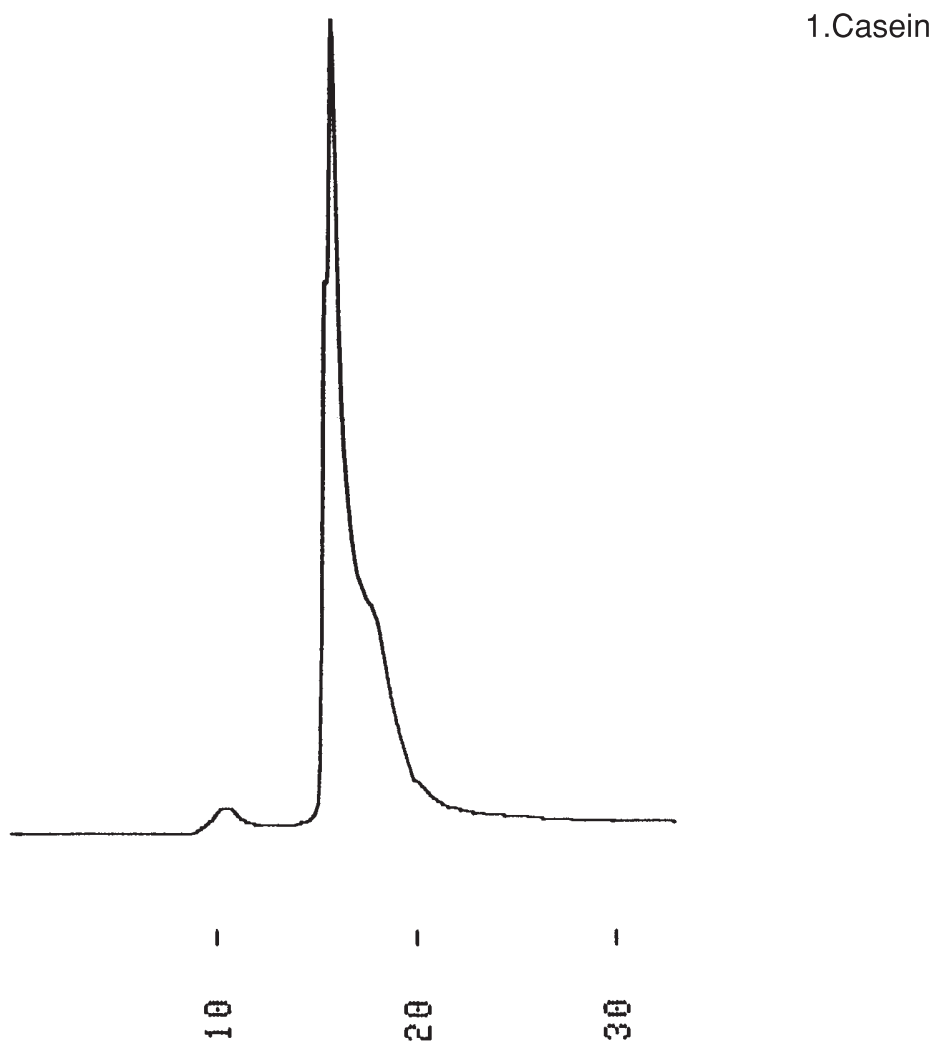

Column	: YMC-Pack Diol-300 500 × 8.0 mm I.D.
Eluent	: 0.1M KH <sub>2</sub> PO <sub>4</sub> -K <sub>2</sub> HPO <sub>4</sub> (pH 7.0) containing 0.2M NaCl/acetonitrile (90/10)
Flow rate	: 1.0 mL/min
Temperature	: ambient (22 °C)
Detectio	: UV at 215 nm, 0.16 AUFS
Injection	: 15 μL (1.0 mg/mL)