


酸化型および還元型グルタチオン
Glutathione, oxidized and reduced form

H900710D

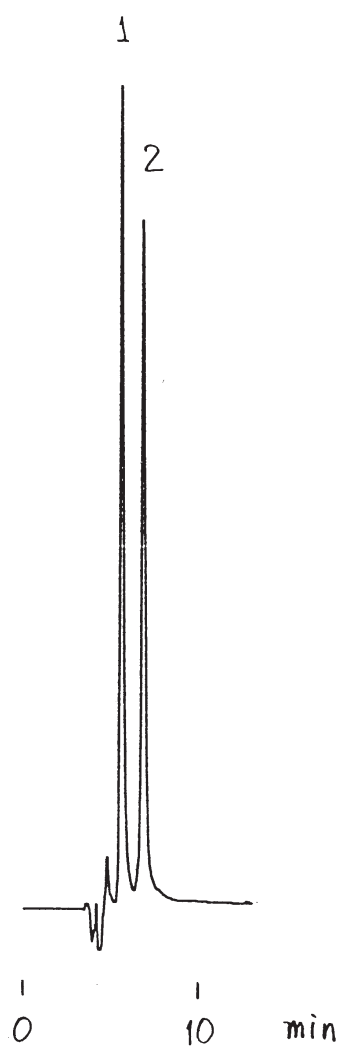

1. Glutathione (GSH)
2. Glutathione disulfide (GSSG)

Column	: YMC-Pack ODS-A (5 μ m, 12nm) 150 \times 4.6 mm I.D.
Eluent	: 0.1M sodium perchlorate in A A) acetonitrile/water/TFA (5/95/0.05)
Flow rate	: 0.5 mL/min
Temperature	: 30 $^{\circ}$ C
Detection	: UV at 220 nm, 0.64 AUFS
Injection	: 10 μ L