
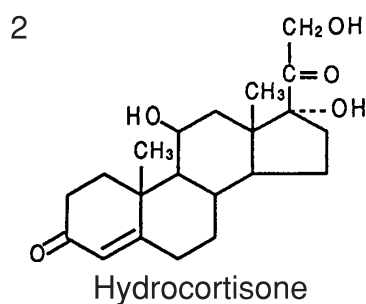
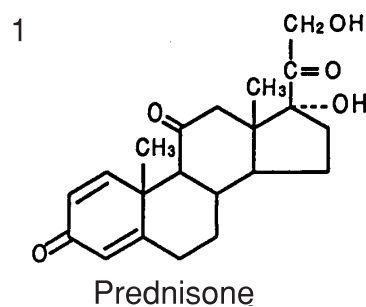


ヒドロコルチゾンとプレドニゾン  
Hydrocortisone and prednisone

K940530C

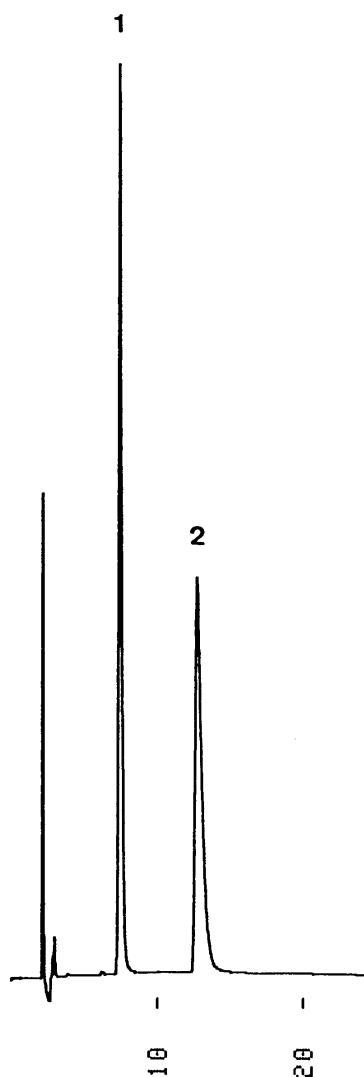

Column : YMC-Pack SIL (5  $\mu$  m, 12nm)  
150  $\times$  4.6 mm I.D.  
Eluent : chloroform/2-propanol/acetic acid (95/5/1)  
Flow rate : 1.0 mL/min  
Temperature : 30  $^{\circ}$ C  
Detection : UV at 250 nm, 0.32 AUFS  
Injection : 10  $\mu$ L (0.25 mg/mL)