


タンパク質
Proteins

M180531BB

- | | |
|-----------------------|--------------|
| 1. Thyroglobulin | (MW 670,000) |
| 2. γ -Globulin | (MW 158,000) |
| 3. Ovalbumin | (MW 45,000) |
| 4. Myoglobin | (MW 17,000) |
| 5. Vitamin B12 | (MW 1,355) |

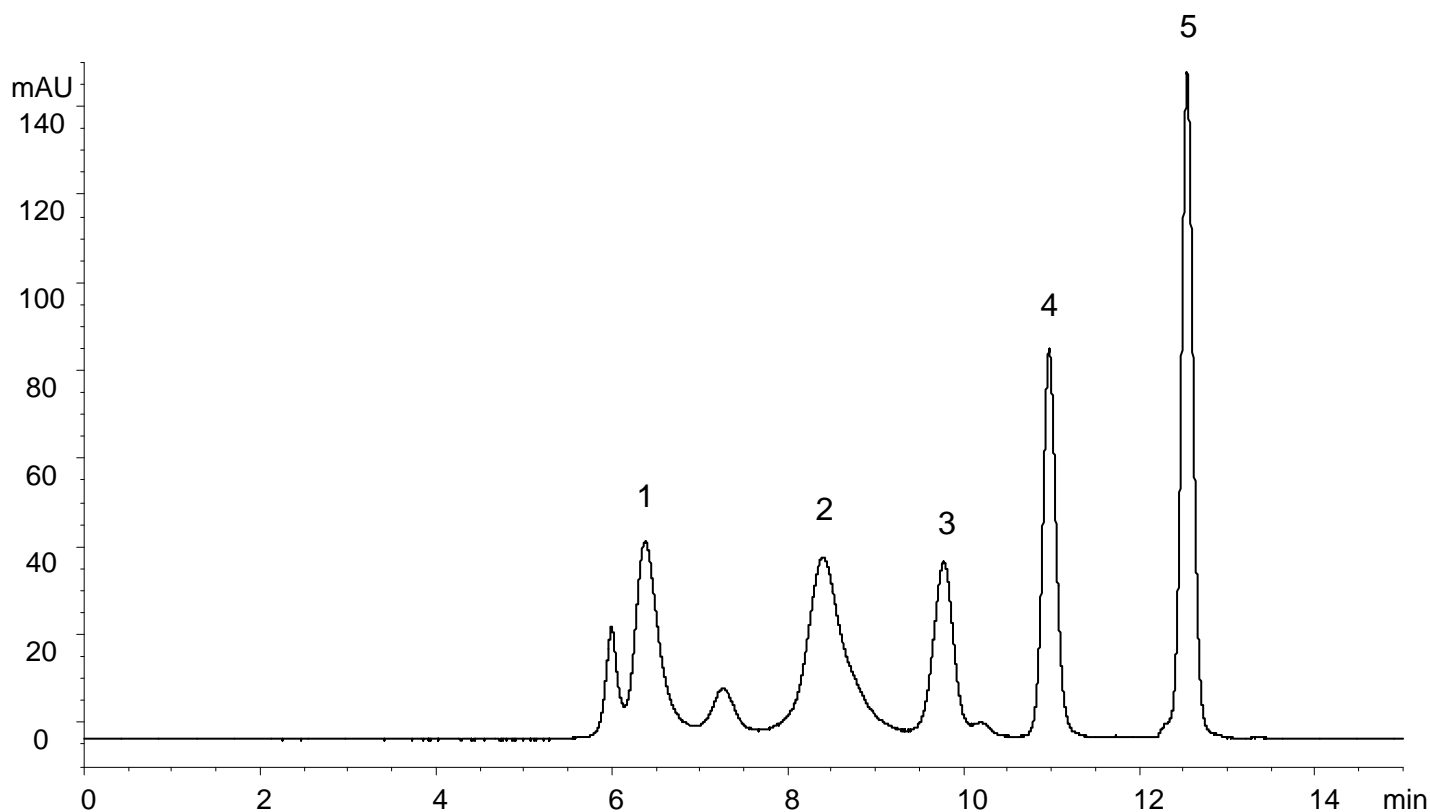

Column : YMC-SEC MAB (3 μ m, 25 nm)
300 X 4.6 mm I.D.

Eluent : 0.1 M NaH₂PO₄-Na₂HPO₄ (pH 7.0) containing 0.2 M NaCl

Flow rate : 0.33 mL/min

Temperature : 25°C

Detection : UV at 280 nm

Injection : 10 μ L (0.14~1.30 mg/mL)