


マルトオリゴ糖(G11~G15)  
Maltooligosaccharides(G11~G15)

M910604K

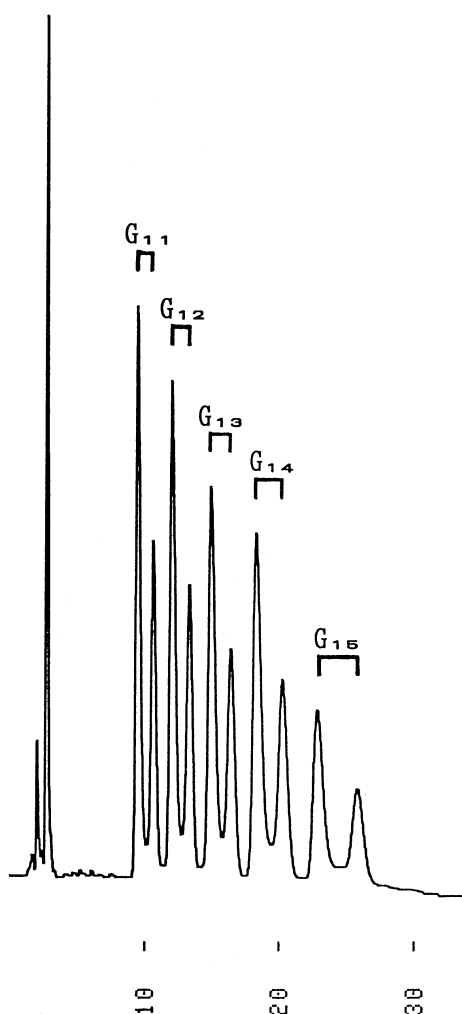

1. Maltoundecaose (G11)
2. Maltododecaose (G12)
3. Maltotridecaose (G13)
4. Maltotetradecaose (G14)
5. Maltopentadecaose (G15)

Column	: YMC-Pack ODS-A (5 $\mu$ m, 30nm) 150 $\times$ 4.6 mm I.D.
Eluent	: water
Flow rate	: 1.0 mL/min
Temperature	: 20 $^{\circ}$ C
Detection	: RI, $8 \times 10^{-6}$ RIU/FS
Injection	: 25 $\mu$ L (1 mg/mL)
Sample	: Manufactured by Nakano Vinegar Co., Ltd., Nakano Central Research Institute