

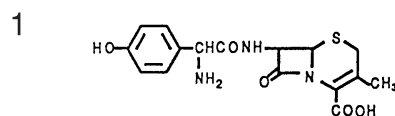


セフェム系抗生物質  
Cephalosporin antibiotics


O940927G

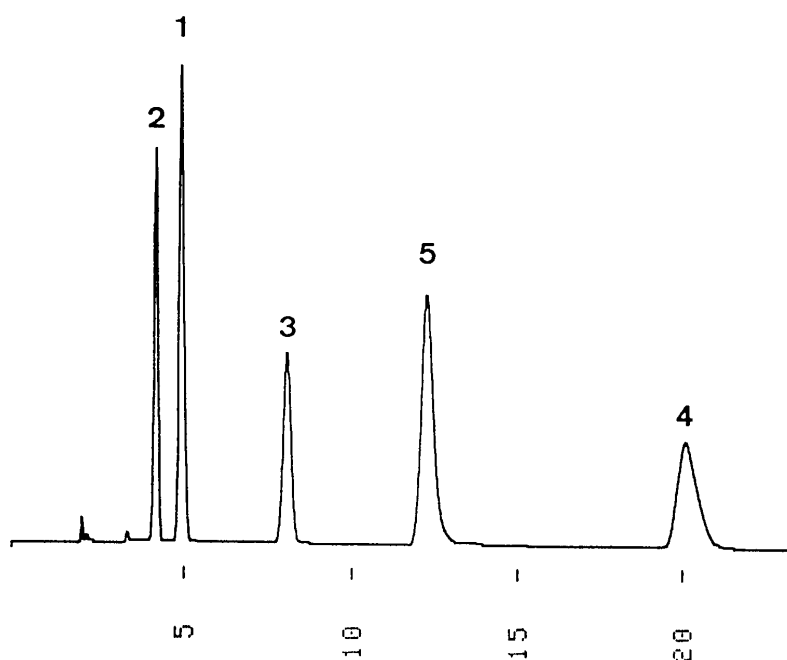


Cefadroxil

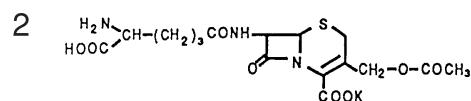


Cephalosporin C potassium

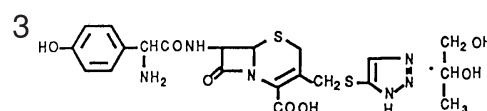

Cefatrizine propylene glycol

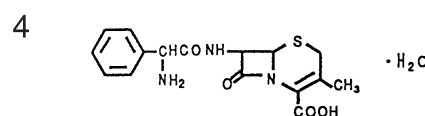

Cephalexin hydrate

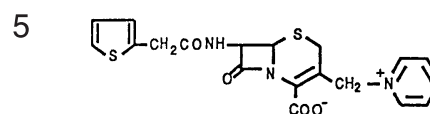

Cephaloridine

Column : YMC-Pack ODS-AQ (5  $\mu$ m, 12nm)  
150  $\times$  4.6 mm I.D.  
Eluent : methanol/water/acetic acid (10/85/5)  
Flow rate : 1.0 mL/min  
Temperature : 37  $^{\circ}$ C  
Detection : UV at 260 nm, 0.16 AUFS  
Injection : 6  $\mu$ L (0.2 mg/mL)