
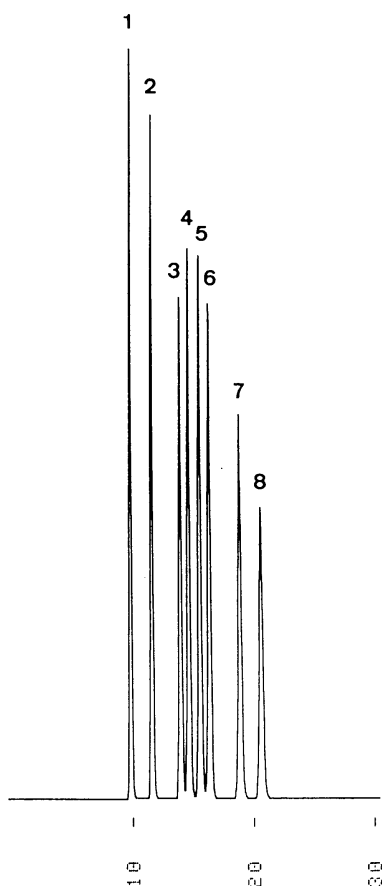
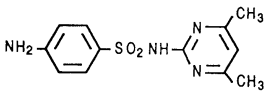
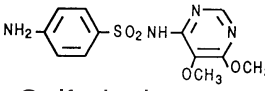
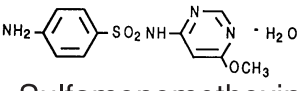
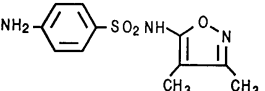
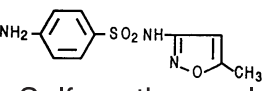
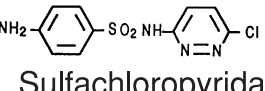
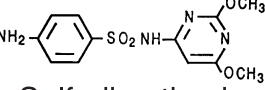


サルファ剤
Sulfonamides

O950302C

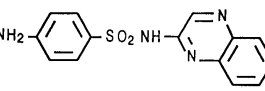

- 1 
Sulfamethazine
- 2 
Sulfadoxine
- 3 
Sulfamonomethoxine hydrate
- 4 
Sulfisoxazole
- 5 
Sulfamethoxazole
- 6 
Sulfachloropyridazine
- 7 
Sulfadimethoxine
- 8 
Sulfaquinoxaline sodium

Column : YMC-Pack ODS-A (5 μ m, 12nm)
250 \times 4.6 mm I.D.
Eluent : THF/water/acetic acid (30/70/1)
Flow rate : 0.5 mL/min
Temperature : 37 $^{\circ}$ C
Detection : UV at 275 nm, 0.16 AUFS
Injection : 5 μ L (42 μ g/mL)