


siRNA duplex under non-denaturing condition

P221202D03

siRNA

5'-CGU ACG CGG AAU ACU UCG AdTdT-3'  
3'-dTdTGCA UGC GCC UUA UGA AGC U-5'

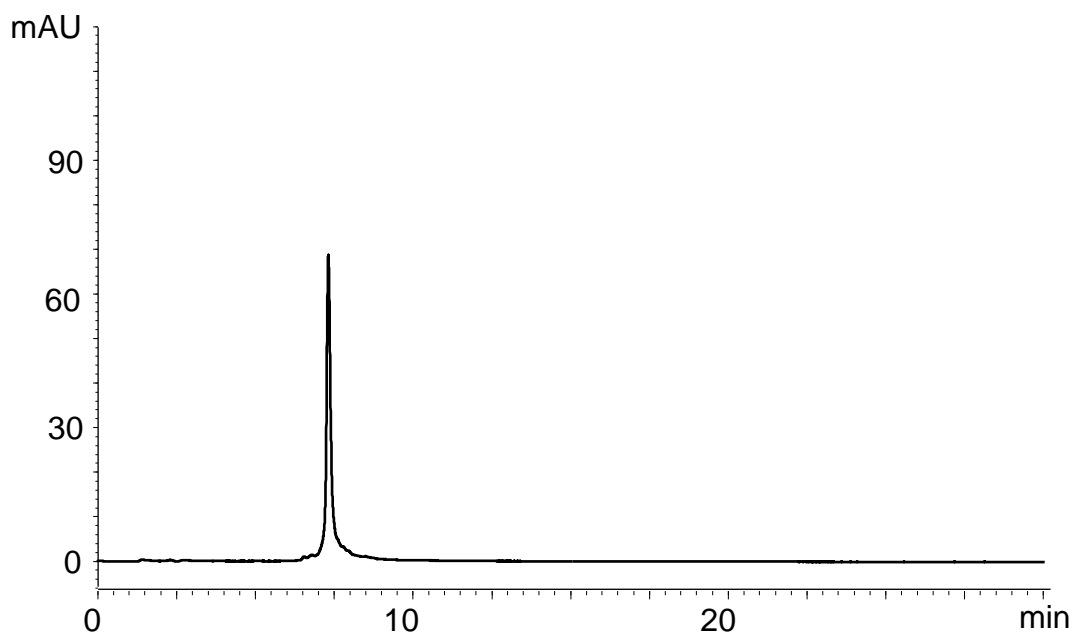

Column	: BioPro IEX QF (5 $\mu$ m) 100 X 4.6 mmI.D.
Eluent	: A) 20 mM Tris-HCl (pH 8.1) B) 20 mM Tris-HCl (pH 8.1) containing 1.0 M NaClO <sub>4</sub> 25-40%B (0-30 min)
Flow rate	: 0.5 mL/min
Temperature	: 25°C
Detection	: UV at 260 nm
Injection	: 4 $\mu$ L (5 nmol/mL)