


有機酸  
Organic acids

Q091022D

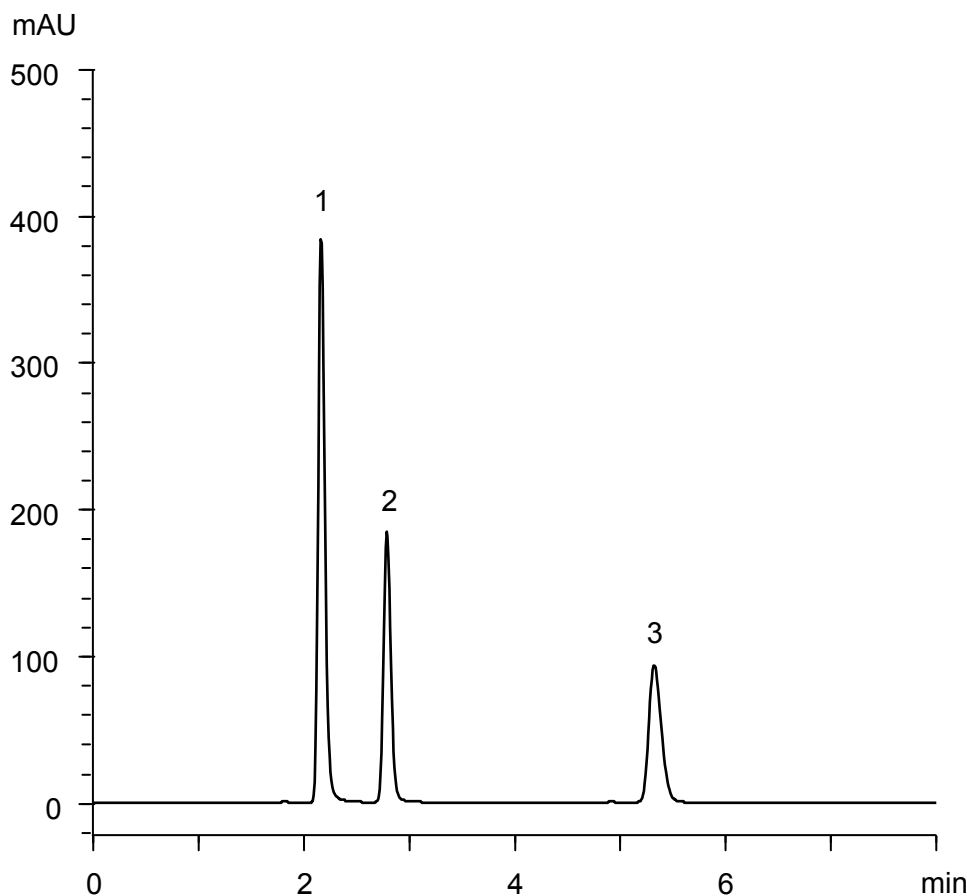

1. Formic acid
2. Acetic acid
3. Propionic acid

Column : YMC-Triart C18 (5  $\mu$ m, 12 nm)  
150 X 3.0 mm I.D.  
Eluent : acetonitrile/0.1% H<sub>3</sub>PO<sub>4</sub> (5/95)  
Flow rate : 0.425 mL/min  
Temperature : 37°C  
Detection : UV at 210 nm  
Injection : 2  $\mu$ L (5  $\mu$ L/mL)