


ペプチド (MW 556 - 3,465)  
Peptides (MW 556 - 3,465)

R091009C

- |                        |            |
|------------------------|------------|
| 1. Oxytocin            | (MW 1,007) |
| 2. Met-Enkephalin      | (MW 574)   |
| 3. Leu-Enkephalin      | (MW 556)   |
| 4. Neurotensin         | (MW 1,673) |
| 5. $\gamma$ -Endorphin | (MW 1,859) |
| 6. $\beta$ -Endorphin  | (MW 3,465) |

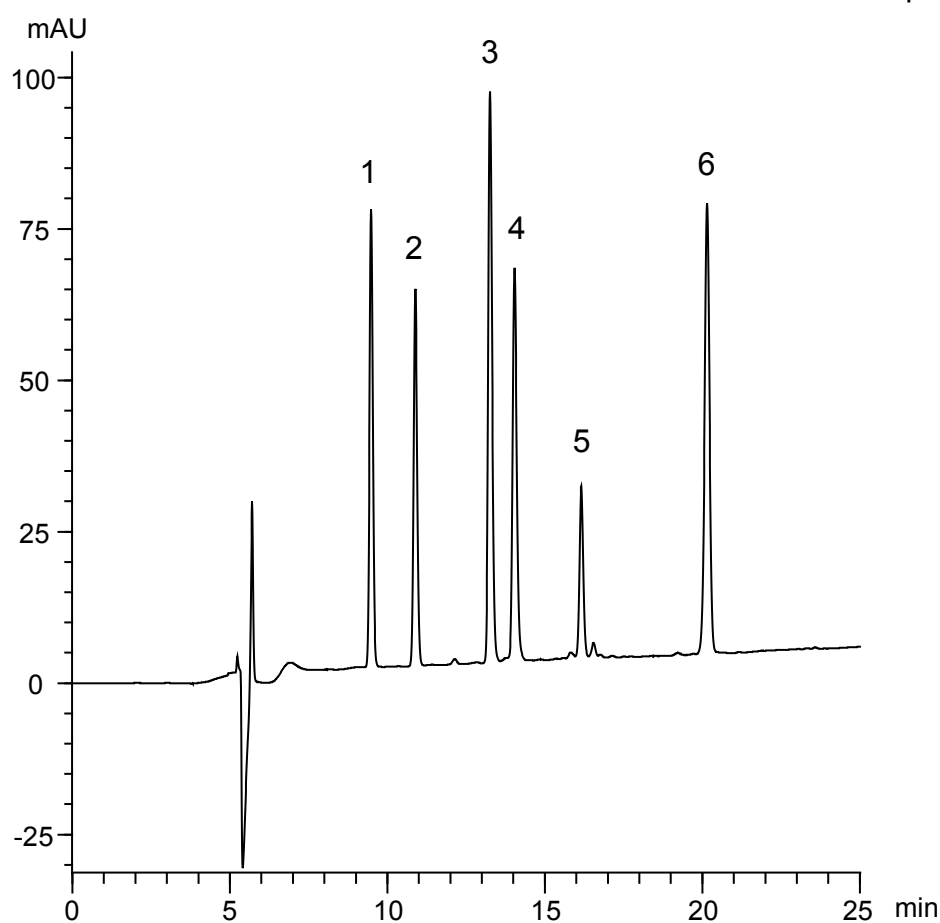

Column : YMC-Triart C18 (5  $\mu$ m, 12 nm)  
150 X 2.0 mm I.D.  
Eluent : A) water/TFA (100/0.1)  
B) acetonitrile/TFA (100/0.1)  
20-45%B (0-25 min)  
Flow rate : 0.2 mL/min  
Temperature : 37°C  
Detection : UV at 220 nm  
Injection : 2  $\mu$ L (0.075~0.25 mg/mL)