



アムロジピンベシル酸塩 (第十五改正日本薬局方第一追補記載条件)  
Amlodipine besilate (The Japanese Pharmacopoeia 15th Edition Supplement I) R100122G

## Standard solution

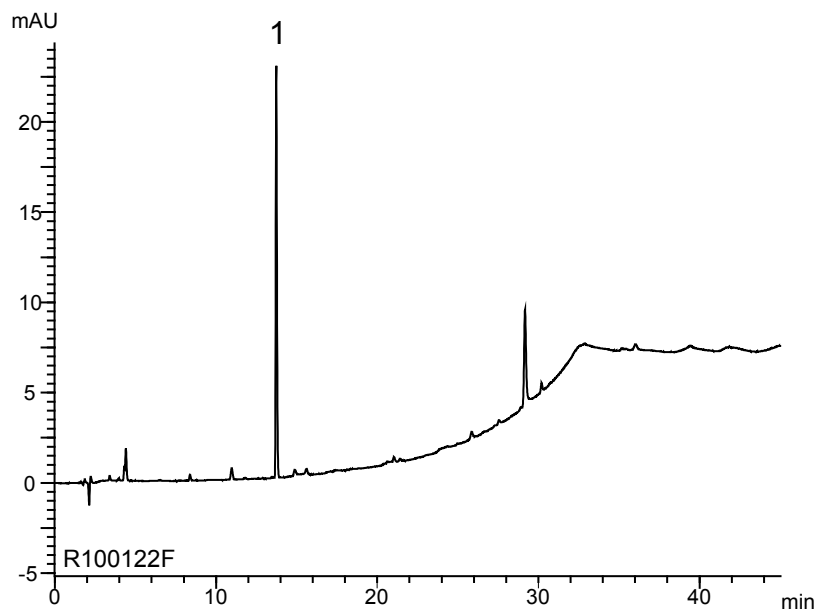


1.

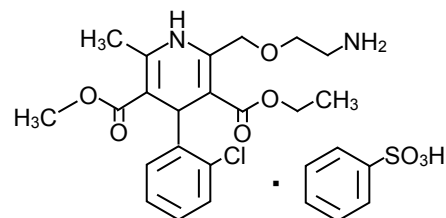


Amlodipine besilate

2.

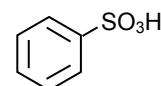

Benzenesulfonic acid

## Sample solution

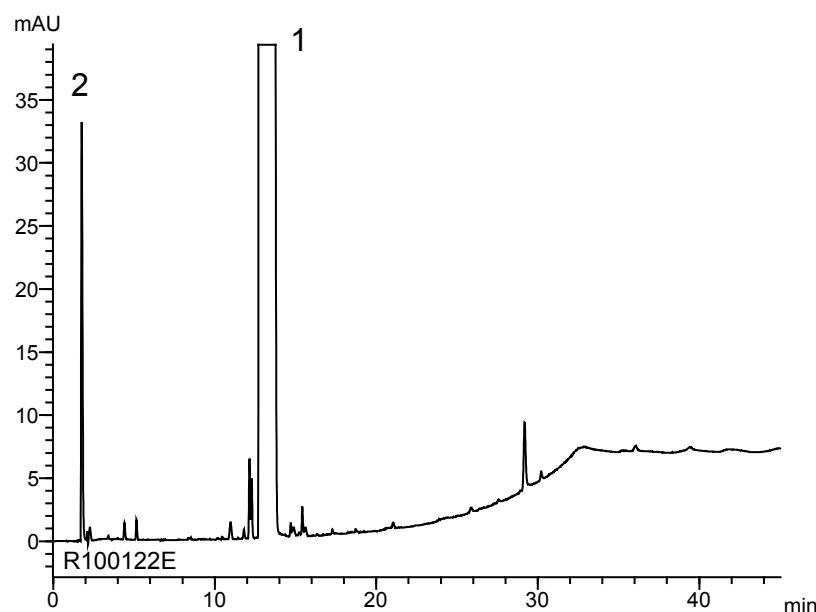

Column : YMC-Triart C18 (3  $\mu$ m, 12 nm)  
150 X 4.6 mm I.D.  
Eluent : A) water/TFA (5000/1)  
B) acetonitrile/TFA (5000/1)  
20-80%B (0-30 min), 80%B (30-45 min)  
Flow rate : 1.0 mL/min  
Temperature : 35°C  
Detection : UV at 237 nm  
Injection : 10  $\mu$ L

(第十五改正日本薬局方第一追補記載条件；純度試験)