
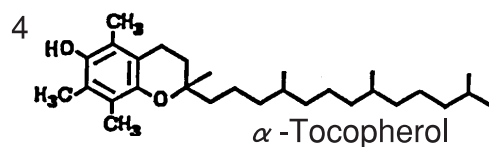
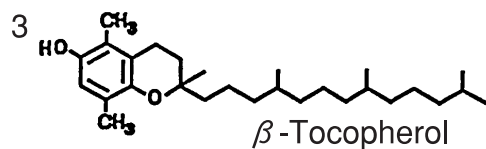
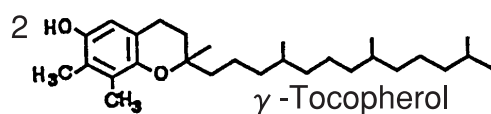
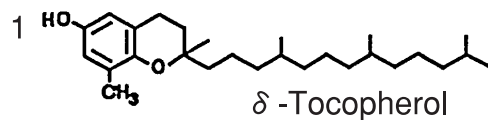


トコフェロール  
Tocopherols

S921216C

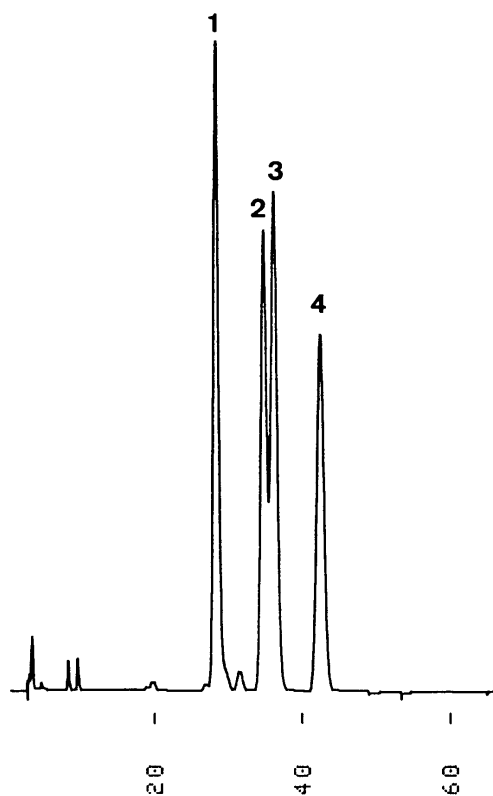

Column : YMC-Pack ODS-A (5  $\mu$ m, 12nm)  
150  $\times$  4.6 mm I.D.  
Eluent : 2-propanol/water (65/35)  
Flow rate : 0.7 mL/min  
Temperature : 30  $^{\circ}$ C  
Detection : UV at 295 nm, 0.04 AUFS  
Injection : 7  $\mu$ L (0.25 mg/mL)  
Ref.) Y.Satomura, M. Kimura and Y. Itokawa, J.Chromatogr., 625,  
372-376 (1992)