

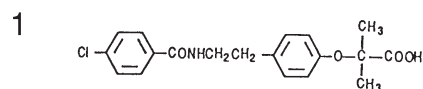


抗高血脂症剂
Hypolipidemic drugs


S930706C

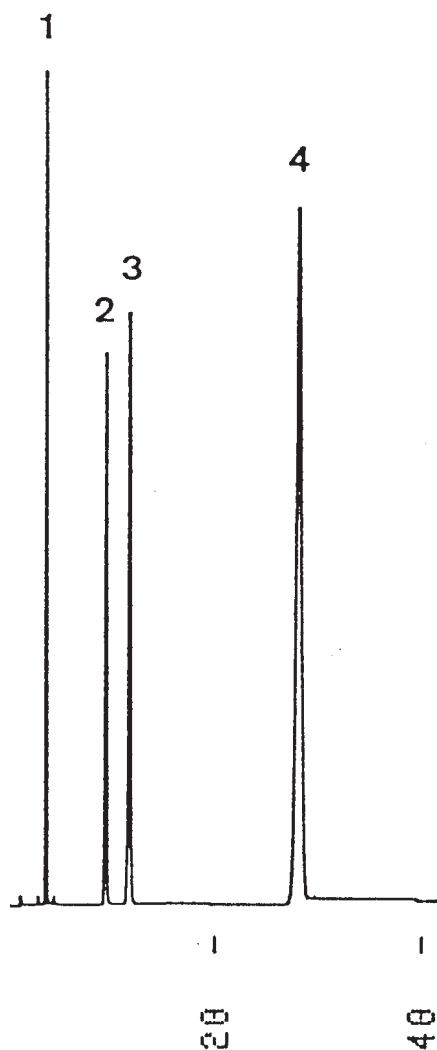


Bezafibrate

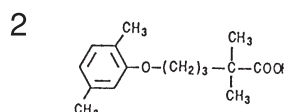

Gemfibrozil

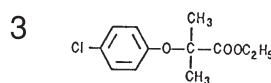

Clofibrate

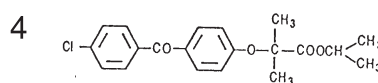

Fenofibrate

Column : J'sphere ODS-M80 (4 μ m, 8nm)
150 \times 4.6 mm I.D.
Eluent : acetonitrile/water/acetic acid (60/40/1)
Flow rate : 1.0 mL/min
Temperature : 30 $^{\circ}$ C
Detection : UV at 275 nm, 0.064 AUFS
Injection : 8 μ L (0.05~0.25 mg/mL)