


水溶性ビタミン
Water-soluble vitamins

S930712E

1. Pyridoxine HCl
2. Nicotinamide
3. Benfotiamine
4. Cyanocobalamin
5. Riboflavin

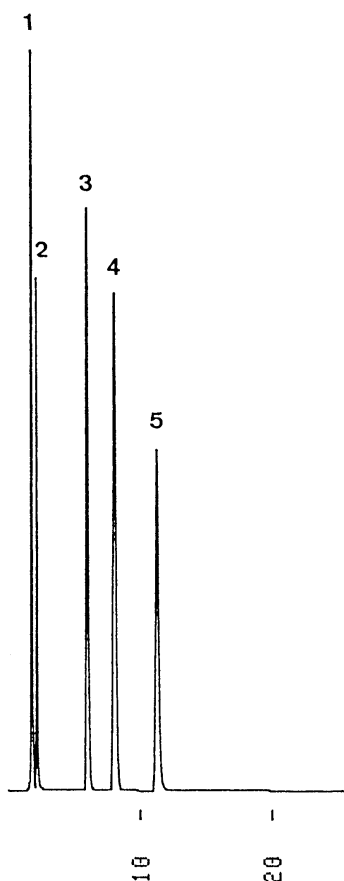

Column : J'sphere ODS-H80 (4 μ m, 8nm)
150 \times 4.6 mm I.D.
Eluent : acetonitrile/50mM (NH₄)₂HPO₄ (10/90)
Flow rate : 1.0 mL/min
Temperature : 30 $^{\circ}$ C
Detection : UV at 254 nm, 0.13 AUFS
Injection : 5 μ L (0.04~0.48 mg/mL)