
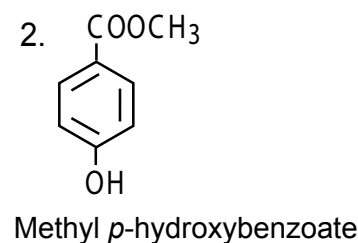
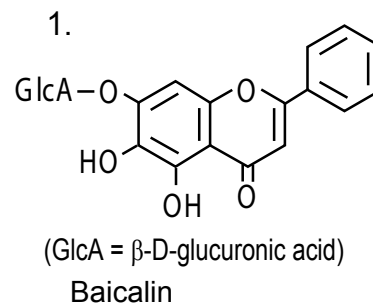


バイカリン  
Baicalin

S961105C

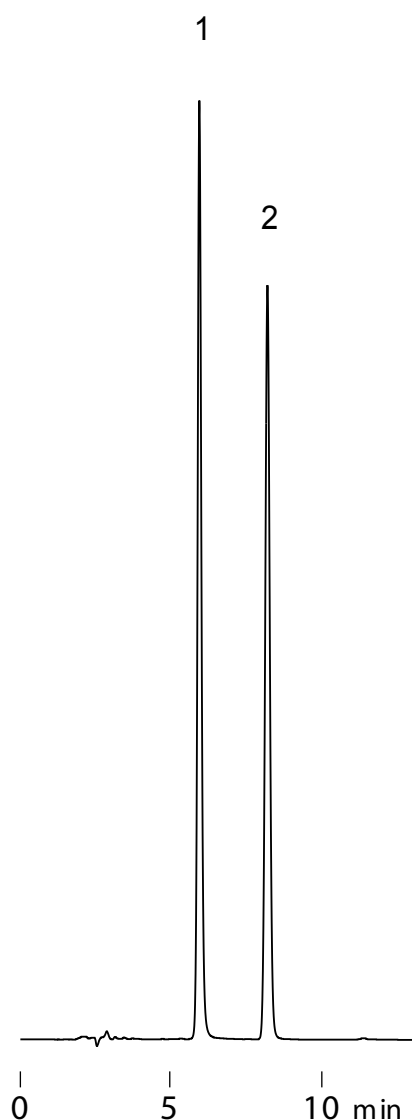

Column : YMC-Pack Pro C18 (5  $\mu$ m, 12 nm)  
150 X 4.6 mm I.D.

Eluent : acetonitrile/10 mM H<sub>3</sub>PO<sub>4</sub> (30/70)

Flow rate : 0.7 mL/min

Temperature : 37°C

Pressure : 5.7 MPa

Detection : UV at 280 nm, 0.064 AUFS

Injection : 10  $\mu$ L (0.01, 0.02 mg/mL)

(日本薬局方第13改正参照)