
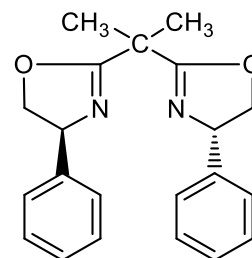



2,2'-イソプロピリデンビス(4-フェニル-2-オキサゾリン)

2,2'-Isopropylidenebis(4-phenyl-2-oxazoline)

T180515F

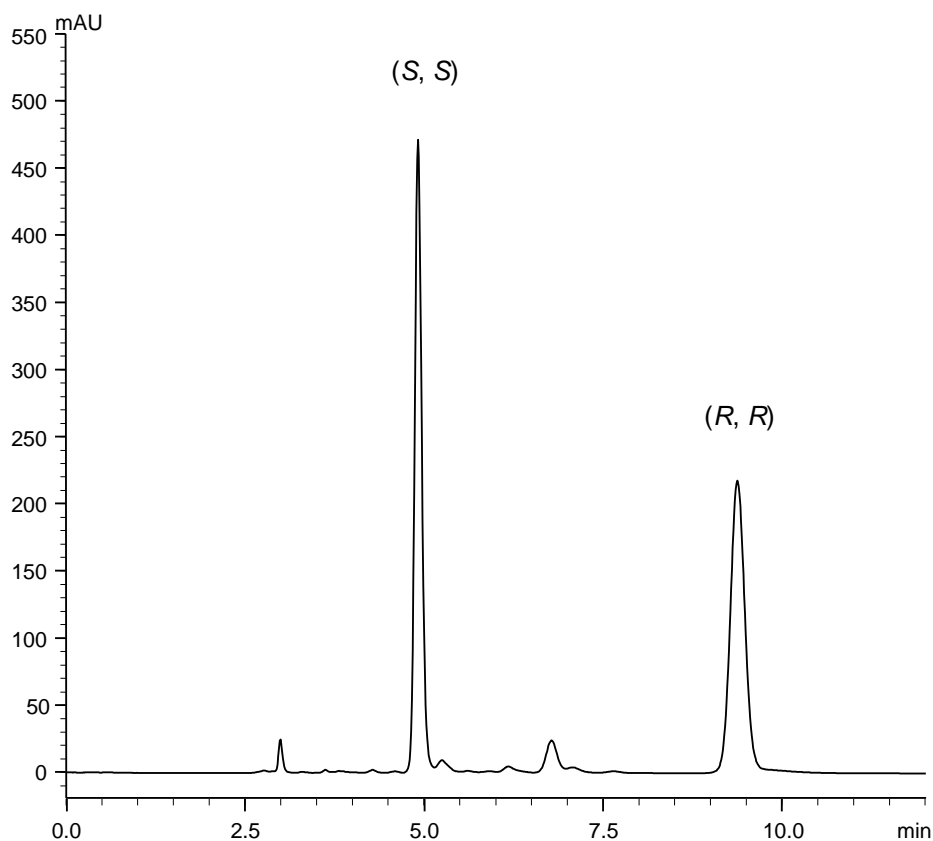

(S, S)

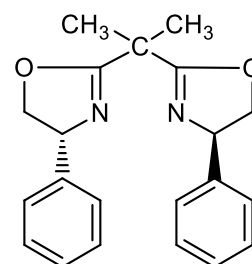

(R, R)

α	3.1
R_s	16.3

Column : CHIRAL ART Amylose-C Neo (5 μ m)
250 X 4.6 mm I.D.
Eluent : *n*-hexane/2-propanol (70/30)
Flow rate : 1.0 mL/min
Temperature : 25°C
Detection : UV at 210 nm
Injection : 10 μ L (0.1 mg/mL)