

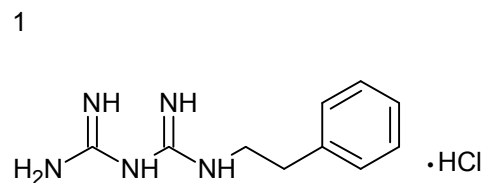


フェンホルミンおよびメトホルミン  
Phenformin and metformin

U120604B

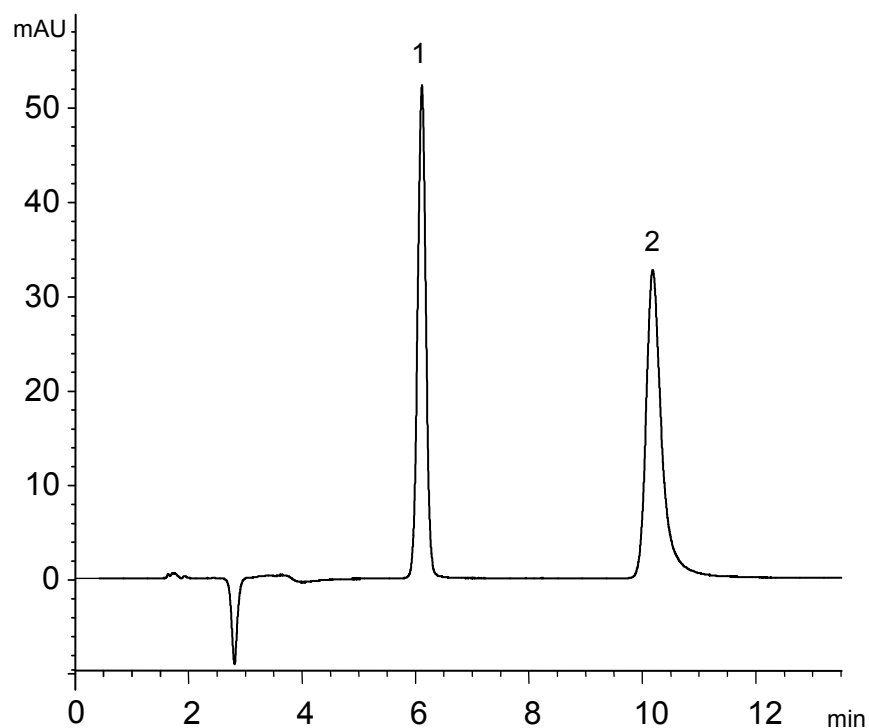

Phenformin hydrochloride

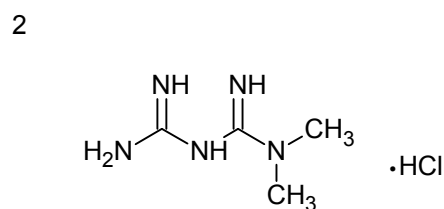

Metformin hydrochloride

Column	: YMC-Triart Diol-HILIC (5 μm, 12 nm) 150 X 3.0 mm I.D.
Eluent	: 100 mM HCOOH-HCOONH <sub>4</sub> (pH 3.7)/acetonitrile (10/90)
Flow rate	: 0.425 mL/min
Temperature	: 25°C
Detection	: UV at 235 nm
Injection	: 5 μL (10 μg/mL)