


合成オリゴヌクレオチド（一本鎖DNA）
Synthetic oligonucleotides (Single-strand DNA)

U190801A

1. 5'-TCATCACACTGAATACCAAT-3' (DNA 20 mer)
2. 5'-**G**TCATCACACTGAATACCAAT-3' (DNA 21 mer)
3. 5'-**A**TCATCACACTGAATACCAAT-3' (DNA 21 mer)
4. 5'-**C**TCATCACACTGAATACCAAT-3' (DNA 21 mer)
5. 5'-**T**TCATCACACTGAATACCAAT-3' (DNA 21 mer)

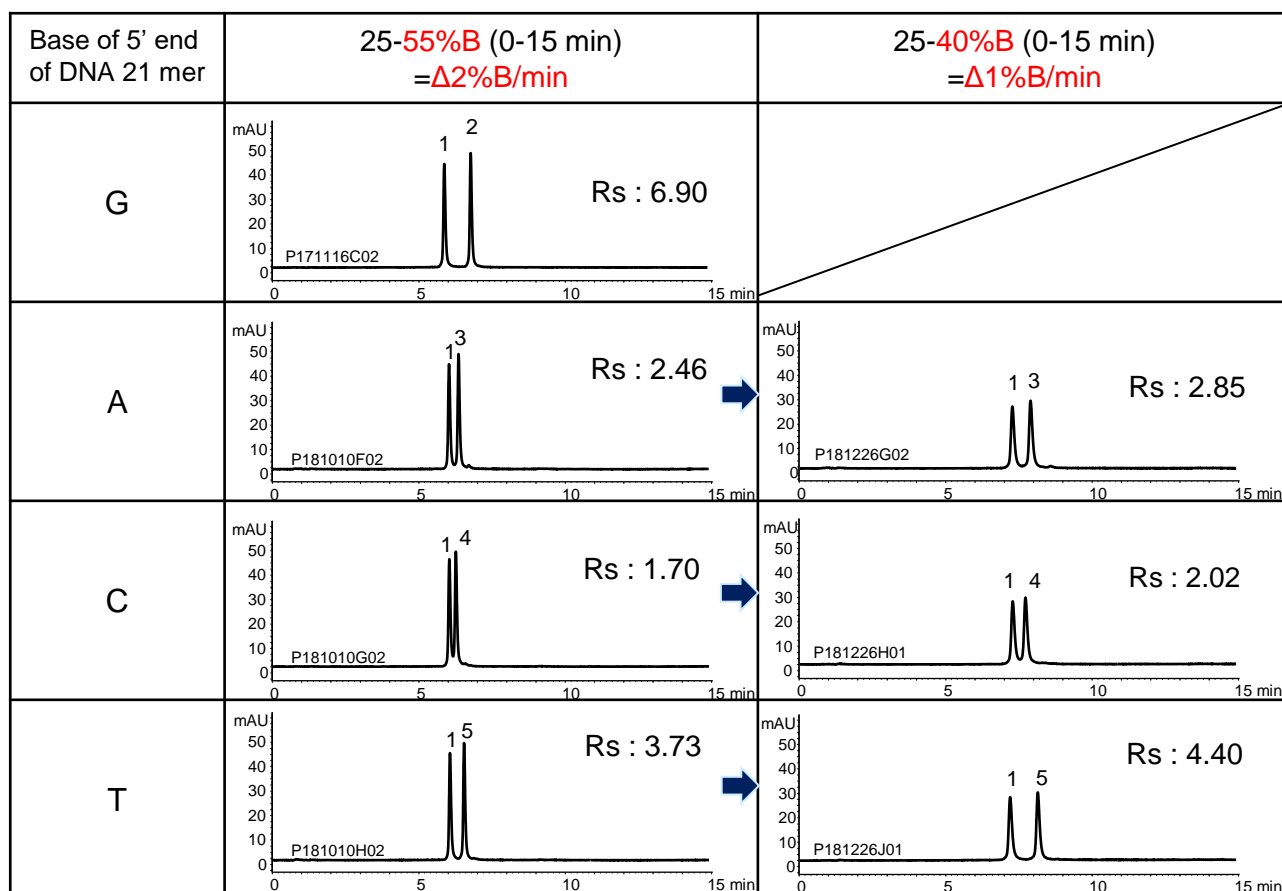

Column : BioPro IEX QF (5 μm)
100 X 4.6 mmI.D.

Eluent : A) 10 mM NaOH
B) 10 mM NaOH containing 1.0 M NaClO₄

Flow rate : 1.0 mL/min

Temperature : 25°C

Detection : UV at 260 nm

Injection : 4 μL (each 5 nmol/mL)