



抗体医薬品
Antibody drugs

U240403A

Bevacizumab

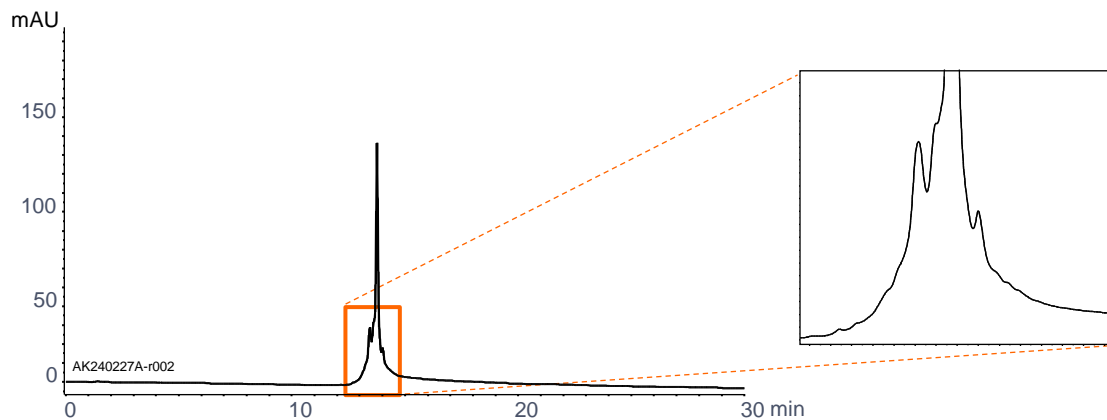

Trastuzumab

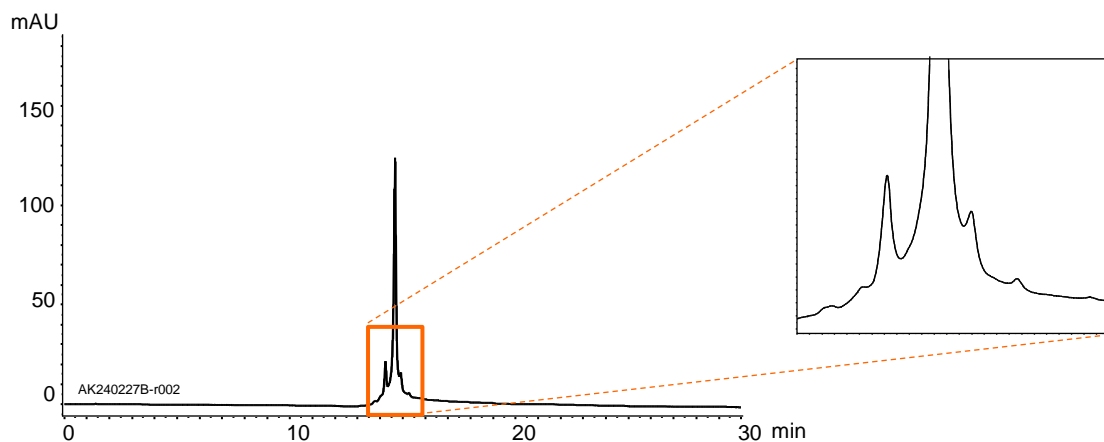

Column : Accura BioPro IEX SF (3 μ m)
100 X 2.1 mmI.D.
Eluent : A) 10 mM MES*-NaOH (pH 6.6)
B) 10 mM MES*-NaOH (pH 6.6) containing 1.0 M NaCl
0-20%B (0-15 min)
Flow rate : 0.1 mL/min
Temperature : 25°C
Detection : UV at 280 nm
Injection : 2 μ L (1 mg/mL)

*2-(*N*-morpholino) ethanesulfonic acid