


トランススチルベンオキサイド
trans-Stilbene oxide


V130812D

1

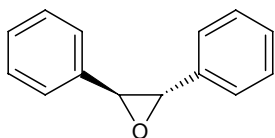


(-)-*S,S*-*trans*-Stilbene oxide

2

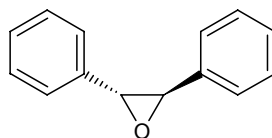

(+)-*R,R*-*trans*-Stilbene oxide

α	2.3
R_s	16.8

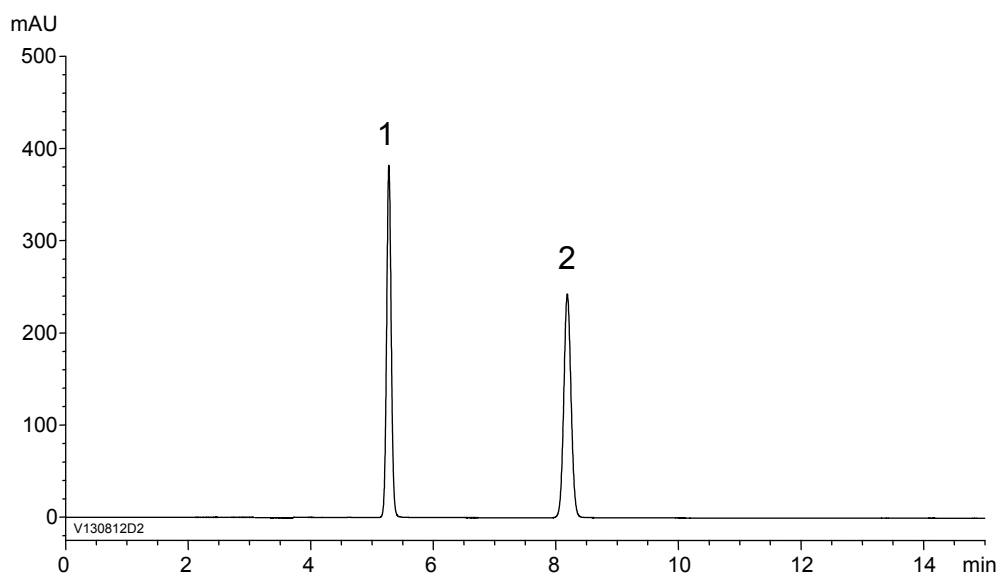

Column : CHIRAL ART Cellulose-C (5 μ m)
250 X 4.6 mm I.D.
Eluent : *n*-hexane/2-propanol (90/10)
Flow rate : 1.0 mL/min
Temperature : 25°C
Detection : UV at 230 nm
Injection : 5 μ L (0.1 mg/mL)