



*N*-CBZ-DL-フェニルアラニン

*N*-CBZ-DL-Phenylalanine

V150403U

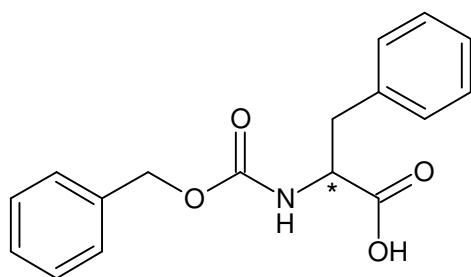

$\alpha$	1.4
$R_s$	4.4

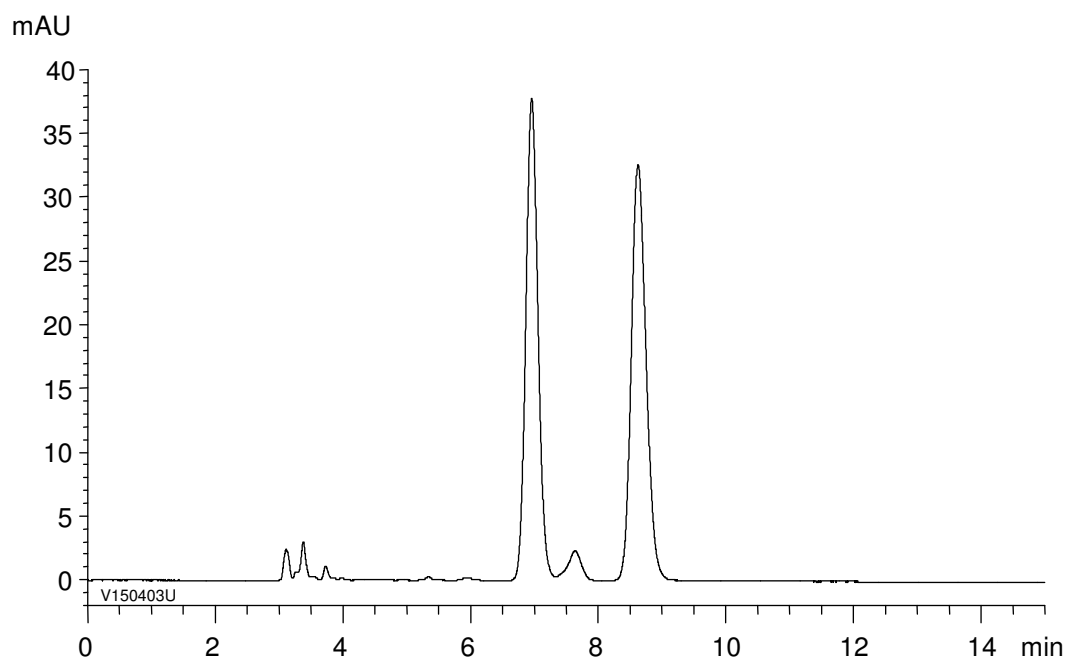

Column : CHIRAL ART Cellulose-SC (5  $\mu$ m)  
250 X 4.6 mm I.D.

Eluent : *n*-hexane/2-propanol/TFA (80/20/0.1)

Flow rate : 1.0 mL/min

Temperature : 25°C

Detection : UV at 254 nm

Injection : 10  $\mu$ L (1 mg/mL)