




モノクローナル抗体フラグメント
Monoclonal antibody fragments

V200204B

Mobile phase containing 0.1 % TFA

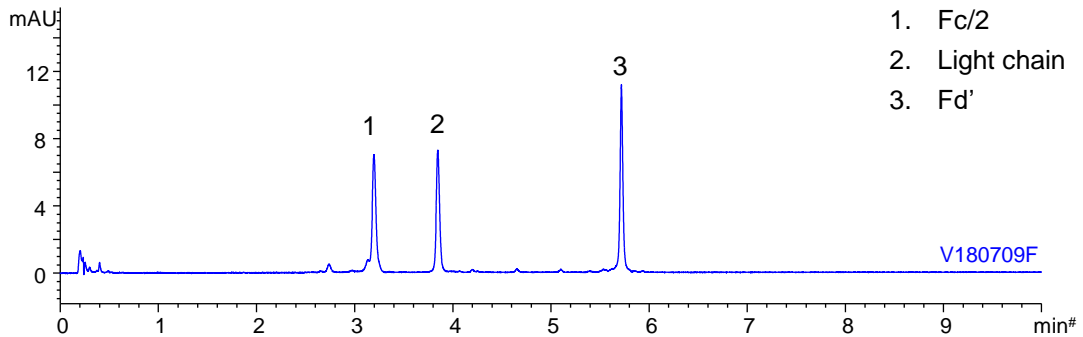

Mobile phase containing 0.1 % formic acid

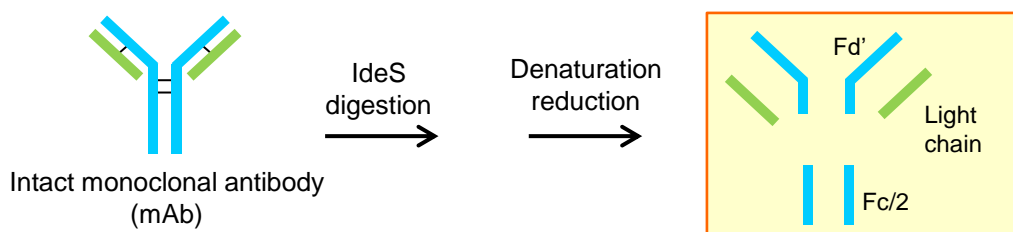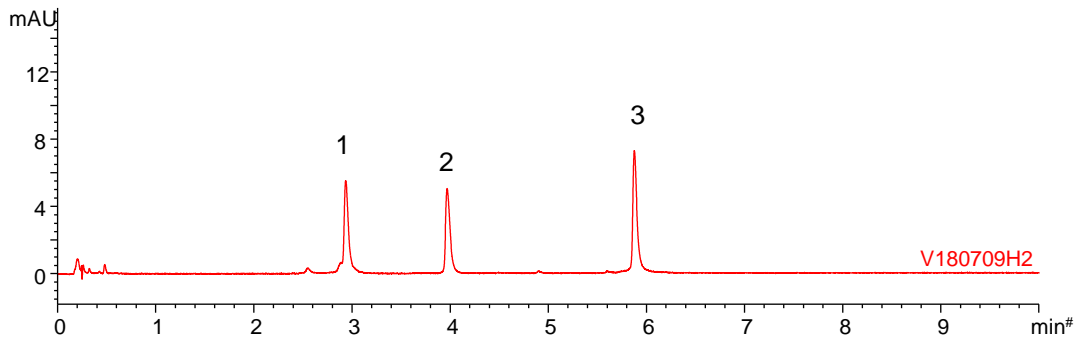

Column	: YMC-Triart Bio C4 (1.9 μ m, 30 nm) 50 X 2.1 mmI.D.
Eluent <TFA>	: A) water/TFA (100/0.1) B) acetonitrile/TFA (100/0.1) 25-50%B (0-10 min), 90%B (10-12.5 min)
Eluent <Formic acid>	: A) water/formic acid (100/0.1) B) acetonitrile/formic acid (100/0.1) 20-45%B (0-10 min), 90%B (10-12.5min)
Flow rate	: 0.4 mL/min
Temperature	: 80°C
Detection	: UV at 280 nm
Injection	: 4 μ L (0.25 mg/mL)
Sample	: Commercial mAb subunit standard