


ペプチドマップ
Peptide mapping

Y890712A

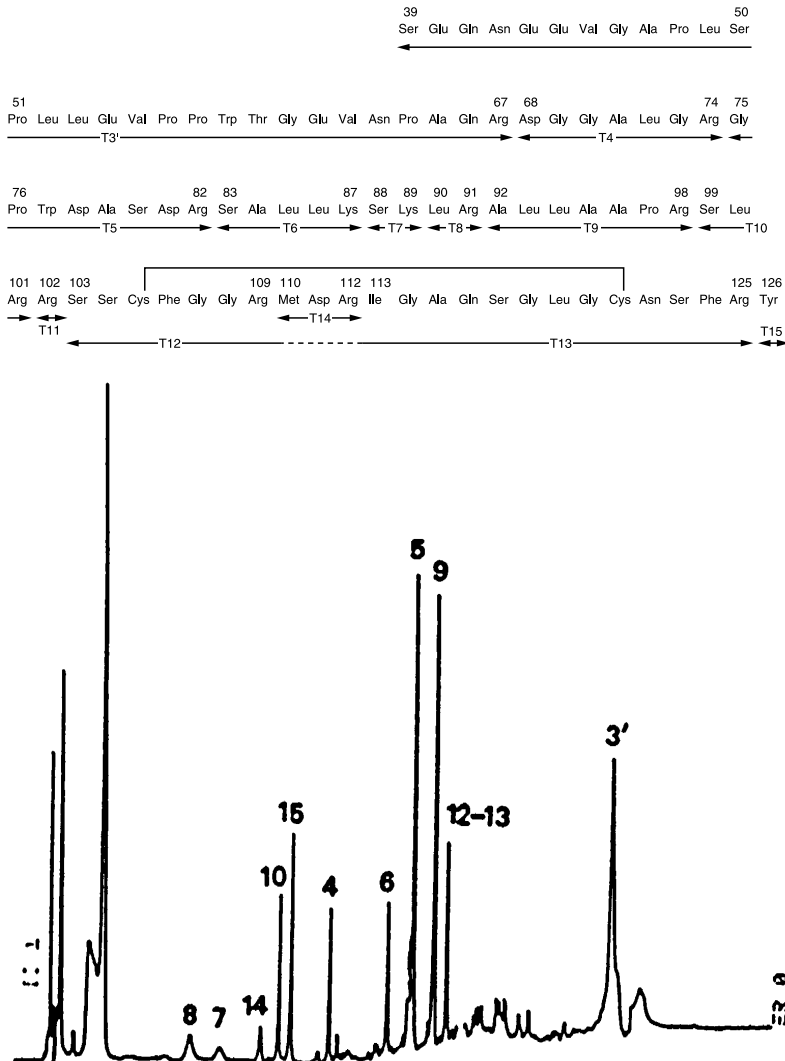

Courtesy of Dr.Nokihara, GBF-Gesellschaft für Biotechnologische Forschung

Column	: YMC-Pack ODS/R (5 μ m, 12nm) 250 \times 4.6 mm I.D.
Eluent	: A)0.01N HCl B)acetonitrile 0%B (0min),0%B (10min),50%B (60min,linear)
Flow rate	: 1.0 mL/min
Detection	: UV 210 nm
Sample	: Tryptic digests of synthetic pCDD 88