
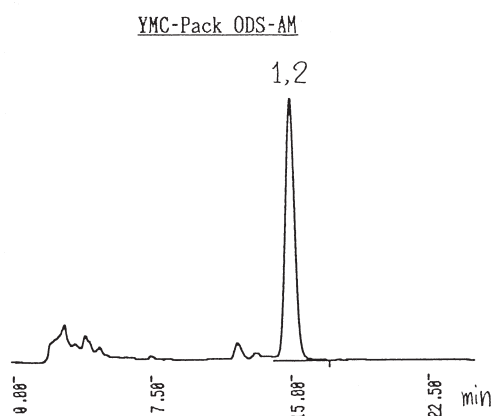
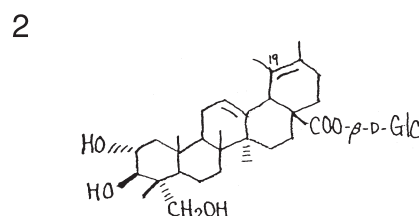
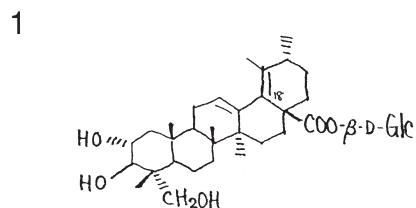


トリテルペン系サポニン

Saponins from roots of rubus suavissimus S. Lee

Y910123

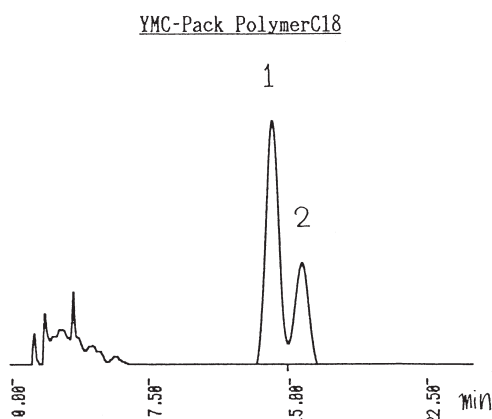

Courtesy of Dr.Ohtani,Hiroshima University,Japan

Column	: YMC-Pack PolymerC18 and YMC-Pack ODS-AM (5 μ m, 12nm) 150 \times 4.6 mm I.D.
Eluent	: methanol/water (45/55) for PolymerC18 methanol/water (50/50) for ODS-AM
Flow rate	: 0.8 mL/min
Temperature	: 50 $^{\circ}$ C
Detection	: UV at 210 nm, 0.32 AUFS
Sample	: 1) β -D-Glucopyranosyl 2 α ,3 β ,23-trihydroxyurs- 12,18-dien-28-oate 2) β -D-Glucopyranosyl 2 α ,3 β ,23-trihydroxyurs- 12,19-dien-28-oate