
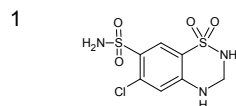



医薬品混合物

Basic drugs


U110801J

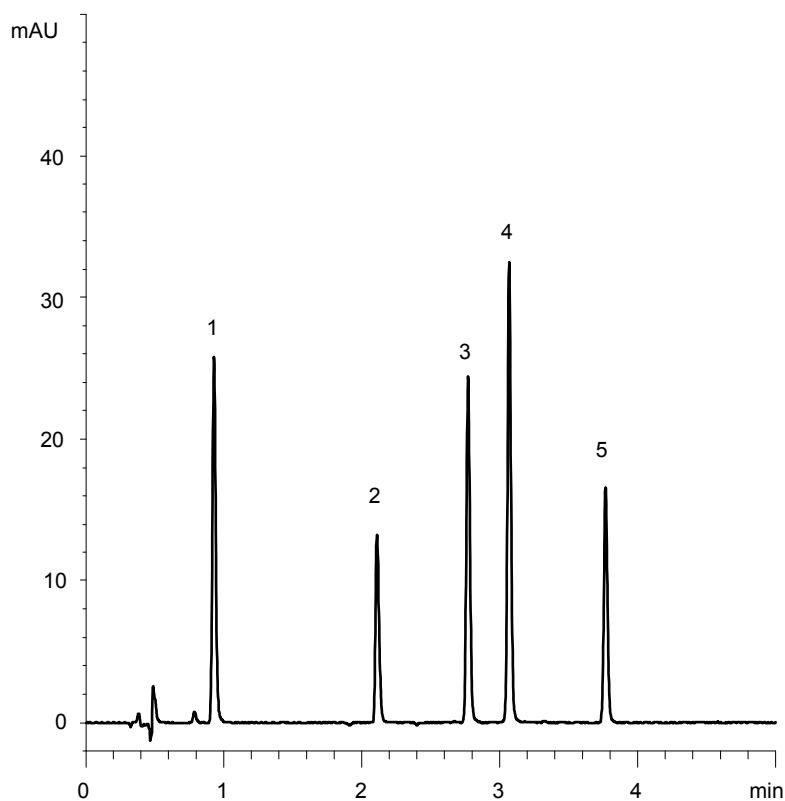


Hydrochlorothiazide

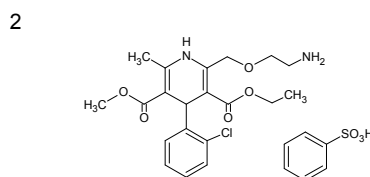


Amlodipine besilate

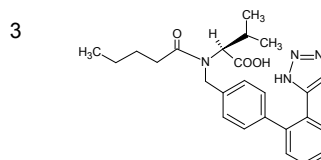

Valsartan

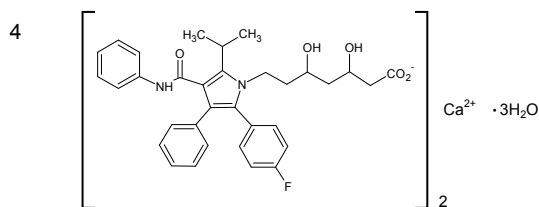

Atorvastatin calcium hydrate

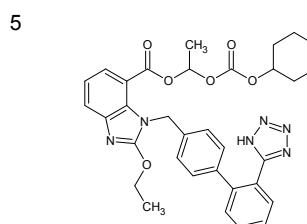

Candesartan cilexetil

Column : YMC-Triart C8 (3 μ m, 12 nm)
50 X 2.0 mmI.D.

Eluent : A) water/formic acid (100/0.1)
B) acetonitrile/formic acid (100/0.1)
10-90%B (0-5 min)

Flow rate : 0.4 mL/min

Temperature : 30°C

Detection : UV at 254 nm

Injection : 2 μ L (10~20 μ g/mL)